



PRODUCT CATALOGUE

MULTIPLE SOLUTIONS FOR
HEALTH CARE PROVIDERS

THERAPIES. HAND IN HAND.

BSNmedical

WOUND CARE

FIXATION - FILMS, TAPES AND BANDAGES

Hypafix® / Cover-Roll® stretch – Non-Woven Adhesive Fixation Sheet	6
Hypafix® transparent – Waterproof Adhesive Transparent Film in Roll	7
Leukomed® T – Waterproof Adhesive Transparent Sterile Film	8
Leukomed® IV – Waterproof Adhesive Transparent Sterile Film for IV	8
Leukoplast® Sleek / Sleek LF – Zinc Oxide Plastic Waterproof Adhesive Tape	9
Elastomull® – Non-Adhesive Fixation Bandage	9
Easifix® – Superior Non-Adhesive Fixation Bandage	10
Easifix® Cohesive – Self-Adhesive Fixation Bandage	10

POST-OP DRESSINGS

Leukomed® – Non-Woven Adhesive Sterile Dressing with Absorbent Pad	11
Leukomed® T plus – Waterproof Adhesive Transparent Sterile Dressing with Absorbent Pad	12
Leukomed® Sorbact® – Antimicrobial Waterproof Transparent Sterile Dressing with Absorbent Pad	13

FIRST AID DRESSINGS

Coverlet® – Stretchable Lightweight Fabric Dressing	13
Covermed® – Non-Woven Adhesive Dressing in Roll Format	14
Coverplast® Classic – Heavyweight Fabric Dressing	14
Coverplast® Barrier – Waterproof & Bacteria Proof Plastic Dressing	15
Coverplast® Standard – Water-Resistant Plastic Dressing	15

IMPREGNATED WOUND CONTACT LAYERS

Cuticell® Classic – Non-Adherent Dressing Impregnated with Paraffin	16
Cuticell® – Non-Adherent Dressing Impregnated with Non-Medicated Ointment	16
Cuticell® Contact – Silicone Wound Contact Layer	17

ANTIMICROBIAL DRESSINGS

Cutimed® Sorbact® – DACC Antimicrobial Dressing	18
---	----

ABSORBENT DRESSINGS

Cutimed® Siltec – Foam Dressing with Silicone Wound Contact Layer	20
Cutimed® Sorbion® – Super-Absorbent Wound Dressing	21
Cutimed® Cavity – Absorbent Foam Dressing for Deep Wounds	22

PRESSURE OFFLOADING

Cutimed® Total Contact Cast Kit – Pressure Offloading Cast	22
--	----

WOUND CLOSURE

Leukosan® Strip – Adhesive Closure Strips	23
---	----

HYDROGEL

Cutimed® Gel – Transparent Amorphous Hydrogel	23
---	----

COMPRESSION THERAPY

BANDAGES

Comprilan® – Short Stretch Compression Bandage	24
JOBST® Compri2 – 2-Layer Bandage System	25
JOBST® Comprifore LF – Multi-Layer Compression Bandage	25
Tensopress® – Long Stretch Compression Bandage	26
Gelocast® – Zinc Paste Compression Bandage (Unna Boot)	26

COMPRESSION SYSTEMS

JOBST® UlcerCARE™ – Ready-to-wear Compression System	27
JOBST® FarrowWrap® – Short Stretch Compression Wrap	28

ACCESSORIES

CompriFoam® – Foam Bandage	30
JOBST® foam – Foam Rubber Pad	30
Artiflex® – Non-Woven Padding Bandage	30

ANTI-EMBOLISM STOCKINGS

Anti-Em/GP™ – Anti-Embolic Stockings with Graduated Compression	31
---	----

ORTHOPEDICS / FRACTURE MANAGEMENT

SYNTHETIC CAST TAPES

Delta-Lite® Plus – Fiberglass Cast Tape	32
Delta-Cast® Conformable – Polyester Cast Tape	33
Delta-Cast® Elite – Polyester Cast Tape	34
Delta-Cast® Prints – Polyester Printed Cast Tape	35
Delta-Cast® Soft – Polyester Semi-Rigid Cast Tape	35

SYNTHETIC SPLINTS

Dynacast® Prelude – Fiberglass Splinting System	36
Delta-Cast® Prelude – Polyester Splinting System	37

PLASTER OF PARIS

Gypsona® – Plaster of Paris Bandages & Slabs	38
Gypsona® Roll – Plaster of Paris Splinting System	38
OCL® – Plaster of Paris Bandages, Slabs & Splinting System	39

UNDERCAST MATERIAL

Delta Terry-Net™ S – Synthetic Terry Cloth Stockinette	39
Delta-Net® / Protouch® – Synthetic Stockinette	40
Protouch® – Cotton Stockinette	40
Delta-Dry® – Water Resistant Cast Padding and Stockinette	41
Specialist® 100 / Specialist® 100 Sterile – Cotton Cast Padding	42
Sof-Rol® / Sof-Rol® Sterile – Natural Cast Padding	42
Protouch® Plus – Synthetic Cast Padding	43
Delta-Rol® – Synthetic Cast Padding	43

ORTHOPEDICS / FRACTURE MANAGEMENT

CAST REMOVAL

American Orthopaedic™ – Cast Removal System	44
Cast Saw Blades	45
Instruments	45

ACCESSORIES

Delta Terry-Net™ – Adhesive Fleece Edger	46
Delta Terry-Net™ – Adhesive Padding	46
Delta Terry-Net™ – Adhesive Hook & Stretch Loop Closures	47
Knee Rest	47
Adjustable Cast Stand	47
Tensoplast® STK – Skin Traction Kit	48
Walking Heels	48
Cast Shoes	48

ORTHOPEDIC SUPPORTS

Actimove Walker Air & Standard – Orthopedic Walker	49
Actimove® Sling – Versatile Upper Limb Support	50

ORTHOPEDICS / PHYSIOTHERAPY - BANDAGES AND TAPES

STRAPPING TAPES

Strappal® – Superior Strapping Tape	51
Leukotape® P – High Adhesive Rigid Strapping Tape	52

ELASTIC ADHESIVE BANDAGES AND TAPES

Tensoplast® / Tensoplast® LF / Tensosport® – Robust Elastic Adhesive Bandage	52
Leukotape® K – Elastic Adhesive Tape for Pain Relief	53
Tensoplast® – Robust Elastic Adhesive Tape	54
Tensoplast® – Athletic Elastic Adhesive Tape	54
Lightplast® Pro – Light Elastic Adhesive Tape	55

COHESIVE BANDAGES

Co-Plus® / Co-Plus® LF – Elastic Cohesive Bandage	55
---	----

ELASTIC NON-ADHESIVE BANDAGES

Econo-San® – Non-Adhesive Elastic Bandage	56
Elastocrepe® – Non-Adhesive Elastic Bandage	56
Elastolex® – Superior Non-Adhesive Elastic Bandage	57

TUBULAR BANDAGES

Tricofix® – Tubular Mesh Bandage	57
Elastogrip® – Tubular Elastic Support Bandage	58
Tensogrip® – Tubular Elastic Support Bandage	58

HYPAFIX® / COVER-ROLL® STRETCH

NON-WOVEN ADHESIVE FIXATION SHEET



Hypafix® / Cover-Roll® stretch are indicated for wide-area fixation and widely used in place of surgical tapes to secure and retain primary dressings on surgical or wound sites. They provide light compression without causing constriction for joints and highly contoured parts of the body.

INDICATIONS

- Fixation of primary dressing on small to large surgical or wound sites
- Ideal use for areas of body requiring flexibility such as elbows and knees
- Ideal base coverage for areas applying the patellofemoral taping technique

FEATURES AND BENEFITS

- Soft, stretchable non-woven polyester adapts well to the body
- Breathable
- Provides light compression without causing constriction
- Skin friendly adhesive holds dressing safely, even in the presence of moisture
- Easily cut to eliminate waste
- Latex Free

Dimensions	Code	Box Contents	Case Contents	Selling Unit
Hypafix®				
2.5 cm x 10 m	7144300	1 roll	-	Box
5 cm x 10 m	7144301	1 roll	-	Box
10 cm x 10 m	7144302	1 roll	-	Box
15 cm x 10 m	7144303	1 roll	-	Box
20 cm x 10 m	7144304	1 roll	-	Box
Cover-Roll® stretch				
5 cm x 1.8 m	4554700	1 roll	12 boxes	Case
10 cm x 1.8 m	4554800	1 roll	12 boxes	Case
15 cm x 1.8 m	4554900	1 roll	12 boxes	Case
20 cm x 1.8 m	4555000	1 roll	12 boxes	Case
30 cm x 1.8 m	4555100	1 roll	12 boxes	Case
5 cm x 9.2 m	4555200	1 roll	12 boxes	Case
10 cm x 9.2 m	4555300	1 roll	12 boxes	Case
15 cm x 9.2 m	4555400	1 roll	12 boxes	Case
20 cm x 9.2 m	4555500	1 roll	12 boxes	Case
30 cm x 9.2 m	4555600	1 roll	12 boxes	Case

HYPAFIX® TRANSPARENT

WATERPROOF ADHESIVE TRANSPARENT FILM IN ROLL



The innovative waterproof film of Hypafix® transparent is very thin and stretchable while maintaining superior strength and durability. It also feels soft and adapts well to the body contours like a second skin.

INDICATIONS

- Skin protection against friction and shearing
- Fixation of dressings, catheters, cannulae, patches and ointment

FEATURES AND BENEFITS

- Easy to apply with the unique red stripe system
- Waterproof and bacteria proof
- Permeable to water vapour and oxygen
- Completely transparent for easy inspection without removing dressing
- Convenient roll format
- Latex Free

Dimensions	Code	Box Contents	Selling Unit
5 cm x 10 m	7237800	1 roll	Box
10 cm x 10 m	7237801	1 roll	Box
15 cm x 10 m	7237802	1 roll	Box



LEUKOMED® T

WATERPROOF ADHESIVE TRANSPARENT STERILE FILM



Leukomed® T is a sterile, breathable, waterproof, transparent film dressing designed to cover catheter or wound sites while providing a barrier to external contaminants such as bacteria, blood and body fluids.

INDICATIONS

- Superficial wounds
- Post-surgical incisions
- Sterile fixation of wound dressings

FEATURES AND BENEFITS

- Easy to apply with the unique red stripe system
- Rounded edges prevent dressing from curling up
- Waterproof and bacteria proof
- Permeable to water vapor and oxygen
- Completely transparent for easy inspection without removing dressing
- Sterile
- Latex Free

Dimensions	Code	Box Contents	Selling Unit
Hospital Pack			
7.2 cm x 5 cm	7238100	50 dressings	Box
8 cm x 10 cm	7238101	50 dressings	Box
10 cm x 12.5 cm	7238108	50 dressings	Box
11 cm x 14 cm	7238102	50 dressings	Box
10 cm x 25 cm	7238114	50 dressings	Box
15 cm x 20 cm	7238109	50 dressings	Box
15 cm x 25 cm	7238110	50 dressings	Box
Retail Pack			
7.2 cm x 5 cm	7238103	5 dressings	Box
8 cm x 10 cm	7238104	5 dressings	Box
10 cm x 12.5 cm	7238111	5 dressings	Box
11 cm x 14 cm	7238105	5 dressings	Box
15 cm x 20 cm	7238112	5 dressings	Box
15 cm x 25 cm	7238113	5 dressings	Box

LEUKOMED® IV

WATERPROOF ADHESIVE TRANSPARENT STERILE FILM FOR IV



When it comes to safety and reliability, choose Leukomed® IV film for securing I.V. catheters and cannulae. This I.V. dressing is completely transparent, bacteria proof and offers excellent adhesive properties.

INDICATIONS

- Fixation of cannulae and peripheral venous catheters
- Puncture sites

FEATURES AND BENEFITS

- Same advantages as Leukomed® T
- V-shaped slit makes fixation easier (except child size)
- Sterile
- Latex Free

Dimensions	Code	Box Contents	Selling Unit
6 cm x 8 cm	7239000	50 dressings	Box
7 cm x 9 cm	7239005	50 dressings	Box
8.5 cm x 11.5 cm (oval)	7239003	50 dressings	Box
4.5 cm x 4.5 cm (child)	7239004	50 dressings	Box



LEUKOPLAST® SLEEK / SLEEK LF

ZINC OXIDE PLASTIC WATERPROOF
ADHESIVE TAPE



Often referred to as pink tape, Leukoplast® Sleek is indicated for any situation where a reliable and conformable waterproof tape is required. Leukoplast® Sleek is easy to unroll, easy to apply and can stay on for hours and days. This tape resists sweat, grease, oil and water.

INDICATIONS

- Fixation of dressings
- Fixation of tubes (E.T., gastric, etc.)

FEATURES AND BENEFITS

- Strong & reliable
- Slight stretch allows the tape to mould to body contours
- Latex Free (Leukoplast® Sleek LF only)

Dimensions	Code	Box Contents	Selling Unit
Leukoplast® Sleek			
2.5 cm x 3 m	7237907	12 rolls	Box
5 cm x 3 m	7237906	6 rolls	Box
Leukoplast® Sleek LF – Latex Free			
2.5 cm x 3 m	7235906	12 rolls	Box
5 cm x 3 m	7235907	6 rolls	Box
2.5 cm x 3 m	7235908	10 spools	Pack
5 cm x 3 m	7235909	5 spools	Pack

ELASTOMULL®

NON-ADHESIVE FIXATION BANDAGE



Elastomull® offers a precision fit, even in areas where there is constant movement such as the elbow or knee. Its superior performance and stretchability helps in controlling bleeding or edema.

INDICATIONS

- Fixation of dressings
- Under padding for compression bandages

FEATURES AND BENEFITS

- Highly elastic
- Particularly skin friendly
- Comfortable for the patient
- Autoclavable
- Virtually lint-free
- Latex Free

Dimensions*	Code	Box Contents	Selling Unit
4 cm x 4 m	2099	20 rolls	Box
6 cm x 4 m	2100	20 rolls	Box
8 cm x 4 m	2101	20 rolls	Box
10 cm x 4 m	2102	20 rolls	Box
12 cm x 4 m	2103	20 rolls	Box

*Stretched

EASIFIX®

SUPERIOR NON-ADHESIVE FIXATION BANDAGE



Easifix® is a durable fixation bandage, highly elastic and very comfortable. Easifix® helps to secure fixation even for awkward physical contours; it causes no restriction to functional movement or blood circulation.

INDICATIONS

- Fixation of dressings
- Under padding for compression bandage

FEATURES AND BENEFITS

- Light, soft and permits air circulation
- Particularly skin friendly
- Autoclavable
- Durable elasticity, high recovery strength
- Latex Free

Dimensions*	Code	Box Contents	Selling Unit
2.5 cm x 4 m	7196101	20 rolls	Bag
5 cm x 4 m	7142910	20 rolls	Bag
7.5 cm x 4 m	7143003	20 rolls	Bag
10 cm x 4 m	7142804	20 rolls	Bag
15 cm x 4 m	7142911	20 rolls	Bag

*Stretched

EASIFIX® COHESIVE

SELF-ADHESIVE FIXATION BANDAGE



Easifix® Cohesive has remarkable body molding properties and conforms to awkward bony prominences without impairing movement. The cohesive coating ensures that the bandage sticks only to itself and not the patient's skin or clothing. The bandage can be reapplied to the same patient, thus permitting routine wound inspection.

INDICATIONS

- Fixation of dressings
- Retention of splinting system
- Retention of compression bandages

FEATURES AND BENEFITS

- Light, soft and permits air circulation
- Economical 20 m roll format
- Lightweight construction
- Non-slip fixation

Dimensions*	Code	Box Contents	Selling Unit
4 cm x 4 m	7143609	10 rolls	Box
6 cm x 4 m	7143610	10 rolls	Box
8 cm x 4 m	7143611	10 rolls	Box
10 cm x 4 m	7143612	10 rolls	Box
12 cm x 4 m	7143613	10 rolls	Box
4 cm x 20 m	7143614	1 roll	Box
6 cm x 20 m	7143615	1 roll	Box
8 cm x 20 m	7143616	1 roll	Box
10 cm x 20 m	7143617	1 roll	Box
12 cm x 20 m	7143618	1 roll	Box

*Stretched

LEUKOMED®

NON-WOVEN ADHESIVE STERILE DRESSING WITH ABSORBENT PAD



Leukomed® adhesive dressing is composed of an acrylic adhesive and non-woven substrate with a highly absorbent wound pad for secreting wounds. The latex-free adhesive enables easy removal, and the non-adherent layer prevents adhesion to the wound.

INDICATIONS

- Post-surgical incisions
- Moderate to heavy drainage wounds
- Acute wounds like cuts, burns and abrasions

FEATURES AND BENEFITS

- Rounded edges prevent dressing from curling up
- Breathable
- Good absorption capacity
- Non-adherent wound pad doesn't stick to wound
- Skin friendly acrylic adhesive
- Highly conformable with two-way stretch
- Sterile
- Latex Free

Dimensions	Code	Box Contents	Selling Unit
Hospital Pack			
7.2 cm x 5 cm	7238000	50 dressings	Box
8 cm x 10 cm	7238001	50 dressings	Box
8 cm x 15 cm	7238002	50 dressings	Box
10 cm x 20 cm	7238003	50 dressings	Box
10 cm x 25 cm	7238004	50 dressings	Box
10 cm x 30 cm	7238005	50 dressings	Box
10 cm x 35 cm	7238006	50 dressings	Box
Retail Pack			
7.2 cm x 5 cm	7238007	5 dressings	Box
8 cm x 10 cm	7238008	5 dressings	Box
8 cm x 15 cm	7238009	5 dressings	Box
10 cm x 20 cm	7238010	5 dressings	Box
10 cm x 25 cm	7238011	5 dressings	Box
10 cm x 30 cm	7238012	5 dressings	Box
10 cm x 35 cm	7238013	5 dressings	Box



LEUKOMED® T PLUS

WATERPROOF ADHESIVE TRANSPARENT STERILE DRESSING WITH ABSORBENT PAD



Leukomed® T plus is a sterile, breathable, waterproof, transparent film dressing with a highly absorbent wound pad for secreting wounds.

INDICATIONS

- Post-surgical incisions
- Moderate to heavy drainage wounds
- Acute wounds like cuts, burns and abrasions

FEATURES AND BENEFITS

- Easy to apply with the unique red stripe system
- Rounded edges prevent dressing from curling up
- Waterproof and bacteria proof
- Permeable to water vapor and oxygen
- Convenient, all-in-one sterile dressing
- Good absorption capacity
- Non-adherent wound pad doesn't stick to wound
- Sterile
- Latex Free

Dimensions	Code	Box Contents	Selling Unit
Hospital Pack			
7.2 cm x 5 cm	7238200	50 dressings	Box
8 cm x 10 cm	7238201	50 dressings	Box
8 cm x 15 cm	7238202	50 dressings	Box
10 cm x 20 cm	7238214	50 dressings	Box
10 cm x 25 cm	7238203	50 dressings	Box
10 cm x 30 cm	7238204	50 dressings	Box
10 cm x 35 cm	7238205	50 dressings	Box
Retail Pack			
7.2 cm x 5 cm	7238206	5 dressings	Box
8 cm x 10 cm	7238207	5 dressings	Box
8 cm x 15 cm	7238208	5 dressings	Box
10 cm x 20 cm	7238215	5 dressings	Box
10 cm x 25 cm	7238209	5 dressings	Box
10 cm x 30 cm	7238210	5 dressings	Box
10 cm x 35 cm	7238211	5 dressings	Box

LEUKOMED® SORBACT®

ANTIMICROBIAL WATERPROOF
TRANSPARENT STERILE DRESSING
WITH ABSORBENT PAD



If surgical circumstances or patient risk factors increase the likelihood of an infection or if the wound is already contaminated, Leukomed® Sorbact® is the right choice. Purely by its physical mode of action, it binds microbes to the pad and helps to prevent critical colonization.

INDICATIONS

- Post-surgical incisions
- Post-operative dehisced wounds and lacerations
- Acute wounds like cuts, burns and abrasions

FEATURES AND BENEFITS

- Same advantages as Leukomed® T plus
- Helps to minimize microbial colonization
- Effectively removes bacteria from the wound
- Sterile
- Latex Free

Dimensions	Code	Box Contents	Selling Unit
5 cm x 7.2 cm	7619900	20 dressings	Box
8 cm x 10 cm	7619901	20 dressings	Box
8 cm x 15 cm	7619902	20 dressings	Box
10 cm x 20 cm	7619903	20 dressings	Box
10 cm x 25 cm	7619904	20 dressings	Box
10 cm x 30 cm	7619905	20 dressings	Box
10 cm x 35 cm	7619906	20 dressings	Box
5 cm x 7.2 cm	7619907	3 dressings	Box
8 cm x 10 cm	7619908	3 dressings	Box
8 cm x 15 cm	7619909	3 dressings	Box

COVERLET®

STRETCHABLE LIGHTWEIGHT
FABRIC DRESSING



Coverlet® fabric adhesive bandages are made to absorb liquid fast - faster than any of the competitors tested. Featuring an extra-large wound pad capable of absorbing 10 times its own weight in water, the 360° adhesive surrounding the island pad helps seal the afflicted area.

INDICATIONS

- Minor acute wounds like cuts, lacerations, abrasions

FEATURES AND BENEFITS

- Available in 11 different shapes and sizes
- Perfect for use on knees, elbows, chins, noses and knuckles
- Elastic fabric moves with the curves and contours of body
- Sterile (except for toe shield and eye occluder)
- Latex Free (except for eye occluder)

	Dimensions	Code	Box Contents (dressings)	Selling Unit
	1.9 cm x 7.6 cm	0023000	100	Box
	2.5 cm x 7.6 cm	0023100	100	Box
	2.2 cm	0030100	100	Box
	2.5 cm x 3.1 cm	0030300	100	Box
	3.8 cm x 5 cm	0033000	100	Box
	5 cm x 7.6 cm	0034000	50	Box
	3.8 cm x 7.6 cm Knuckle	0139000	100	Box
	6.3 cm x 3.5 cm Toe	0080100	100	Box
	3.8 cm x 5.4 cm Fingertip	0130600	100	Box
	5 cm x 6.3 cm Fingertip	0130700	50	Box
	7.6 cm x 7.6 cm	0038500	50	Box
	4.8 cm x 6.3 cm Child eye	4642900	20	6 boxes
	5 cm x 7.5 cm Adult eye	4643000	20	6 boxes

COVERMED®

NON-WOVEN ADHESIVE DRESSING IN ROLL FORMAT



Covermed® is made of highly flexible polyester to allow full freedom of movement, let the skin breathe and dry quickly if wet. Comfortable and safe for the patient, it will stay in place for the duration of treatment.

INDICATIONS

- Minor acute wounds like cuts, lacerations, abrasions
- After blood drawing or vaccination

FEATURES AND BENEFITS

- Breathable
- Good absorption capacity
- Non-adherent wound pad doesn't stick to wound
- Skin-friendly acrylic adhesive
- Highly conformable with two-way stretch
- Latex Free

Dimensions	Code	Box Contents	Selling Unit
4 cm x 5 m	7215100	1 roll	Box
6 cm x 5 m	7215200	1 roll	Box
8 cm x 5 m	7215300	1 roll	Box

COVERPLAST® CLASSIC

HEAVYWEIGHT FABRIC DRESSING



For flexible comfort that stretches with every movement, Coverplast® Classic is a strong and durable dressing that protects and cushions the wound, providing the practical answer to everyday minor injuries.

INDICATIONS

- Minor acute wounds like cuts, lacerations, abrasions
- After blood drawing or vaccination

FEATURES AND BENEFITS

- Lasting and durable adhesive
- Robust substrate made of 100% cotton
- Low adherent absorbent wound pad
- Breathable
- Sterile (except for roll format)

Dimensions	Code	Box Contents	Selling Unit
Roll Format			
4 cm x 5 m	7115104	1 roll	Box
6 cm x 5 m	7115103	1 roll	Box
8 cm x 5 m	7115102	1 roll	Box
Pre-Cut Format			
2.2 cm x 7.2 cm	7260107	100 dressings	Box
2.2 cm x 3.8 cm	7260100	100 dressings	Box
3.8 cm x 3.8 cm	7260101	100 dressings	Box
5 cm x 7.2 cm	7260104	100 dressings	Box
3.8 cm x 7.5 cm Knuckle	7260106	100 dressings	Box
4.5 cm x 5 cm Fingertip	7260120	50 dressings	Box
4.5 cm x 7.5 cm Fingertip	7260105	100 dressings	Box
Doctor's Set	7211806	101 dressings	Box
Bulk			
2.2 cm x 7.2 cm	7260117	10000 dressings	Case
5 cm x 7.2 cm	7260113	5000 dressings	Case
4.5 cm x 5 cm Fingertip	7260115	8000 dressings	Case
4.5 cm x 7.5 cm Fingertip	7260114	5000 dressings	Case
3.8 cm x 7.5 cm Knuckle	7260116	5000 dressings	Case

COVERPLAST® BARRIER

WATERPROOF & BACTERIA PROOF
PLASTIC DRESSING



Coverplast® Barrier is waterproof and provides a barrier to the risk of contamination and secondary bacterial infections. Patients can wash and bathe as normal without compromising the protection.

INDICATIONS

- Minor acute wounds like cuts, lacerations, abrasions
- After blood drawing or vaccination

FEATURES AND BENEFITS

- Waterproof and bacteria proof
- Breathable
- Skin friendly acrylic adhesive
- Sterile
- Latex Free

	Dimensions	Code	Box Contents	Selling Unit
	2.2 cm x 7.2 cm	7214301	100 dressings	Box
	2.2 cm x 3.8 cm	7214304	100 dressings	Box
	2.4 cm	7214303	100 dressings	Box
	3.8 cm x 3.8 cm	7214302	100 dressings	Box
	3.8 cm x 6.3 cm	7214305	100 dressings	Box
	5 cm x 7.2 cm	7214306	100 dressings	Box

COVERPLAST® STANDARD

WATER-RESISTANT PLASTIC
DRESSING



Coverplast® Standard washable adhesive dressing protects minor everyday injuries in wet or dirty environments. The flexible perforated film lets the skin breathe and ensures high wearing comfort.

INDICATIONS

- Minor acute wounds like cuts, lacerations, abrasions
- After blood drawing or vaccination

FEATURES AND BENEFITS

- Flexible and conformable
- Low adherent absorbent wound pad
- Washable
- Water and dirt resistant
- Sterile
- Latex Free

	Dimensions	Code	Box Contents	Pack Contents	Selling Unit
	1.9 cm x 7.2 cm	7259001	100 dressings	6 boxes	Pack

CUTICELL® CLASSIC

NON-ADHERENT DRESSING IMPREGNATED WITH PARAFFIN



Cuticell® Classic is a low-adherent primary wound contact layer dressing made of open weave cotton gauze impregnated with soft paraffin to facilitate dressing change and to support a moist environment for the wound.

INDICATIONS

For superficial exuding wounds such as

- First and second-degree burns
- Lacerations and abrasions
- Donor and recipient skin graft sites

FEATURES AND BENEFITS

- Low-adherent and soothing to minimize pain at dressing change
- Supports a moist environment to help healing
- Sterile (except for dressings not individually wrapped)
- Latex Free

Dimensions	Code	Box Contents	Selling Unit
5 cm x 5 cm	7253801	50 dressings	Box
5 cm x 5 cm	7253800	5 dressings	Box
10 cm x 10 cm	7253803	100 dressings	Box
10 cm x 10 cm	7253802	10 dressings	Box
10 cm x 40 cm	7253804	10 dressings	Box
10 cm x 7 m*	7253806	1 roll	Box
15 cm x 2 m*	7253807	12 rolls	Box

* Not individually wrapped

CUTICELL®

NON-ADHERENT DRESSING IMPREGNATED WITH NON-MEDICATED OINTMENT



Cuticell® is a low-adherent primary wound contact layer made of smooth small meshed acetate fabric impregnated with non-medicated ointment basis. The result: a dressing that ensures the wound is not disturbed and helps maintain a moist environment that supports wound healing.

INDICATIONS

For superficial exuding wounds such as

- First and second-degree burns
- Lacerations and abrasions
- Radiation injuries
- Donor and recipient skin graft sites

FEATURES AND BENEFITS

- Low-adherent and soothing to minimize pain at dressing change
- Supports a moist environment to help healing
- Does not disturb wound healing
- Does not charge the wound with excess ointment
- Simple handling without release paper
- Sterile
- Latex Free

Dimensions	Code	Box Contents	Selling Unit
7.5 cm x 7.5 cm	7253901	50 dressings	Box
7.5 cm x 7.5 cm	7253900	10 dressings	Box
7.5 cm x 20 cm	7253903	50 dressings	Box
7.5 cm x 20 cm	7253902	10 dressings	Box
20 cm x 40 cm	7253904	25 dressings	Box

CUTICELL® CONTACT

SILICONE WOUND CONTACT LAYER



Cuticell® Contact is a silicone wound contact layer that has a perforated transparent structure and is highly conformable. Its optimally perforated structure ensures even highly viscous exudate can easily pass through.

INDICATIONS

For superficial exuding wounds such as

- First and second-degree burns
- Lacerations and abrasions
- Radiation injuries
- Donor and recipient skin graft sites

FEATURES AND BENEFITS

- Low-adherent and soothing to minimize pain at dressing change
- Supports a moist environment to help healing
- Allows exudate to drain freely into an absorbent secondary dressing
- High transparency for easy wound inspection
- Reduced risk of maceration
- Sterile
- Latex Free



Dimensions	Code	Box Contents	Selling Unit
5 cm × 7.5 cm	7268000	5 dressings	Box
7.5 cm × 10 cm	7268001	5 dressings	Box
10 cm × 18 cm	7268002	5 dressings	Box
15 cm × 25 cm	7268003	5 dressings	Box

CUTIMED® SORBACT®

DACC ANTIMICROBIAL DRESSING



Traditional antimicrobial dressings use harmful chemicals to kill bacteria or fungi, leaving behind bacterial cell debris and endotoxins. This can significantly slow or disrupt the healing cycle. Cutimed® Sorbact® uses no chemical agents, only the physical property of hydrophobic interaction, to safely remove bacteria without cell death or negatively impacting the healing process.

INDICATIONS




For all unclean, colonized and infected exuding wounds

- Chronic wounds, such as venous, diabetic or pressure ulcers
- Postoperative and dehisced wounds
- Traumatic wounds
- Wounds after excision of fistulas and abscesses


FEATURES AND BENEFITS

- No healing delays due to bacterial death, cell debris or endotoxin release
- No bacterial resistance
- No cytotoxicity
- No risk of allergies and no contraindications
- Icons on box to choose the right product
- Sterile
- Latex Free

The dressing surface of Cutimed® Sorbact® is coated with DACC, a highly hydrophobic fatty acid derivate which wound pathogens find very attractive. When the dressing is moist, wound bacteria and fungi are naturally drawn to the hydrophobic surface. They become irreversibly bound to the DACC-coated fibers, unable to move back to the wound bed. The bound pathogens are removed with each dressing change, which helps clean contaminated, colonized and infected wounds by reducing the overall bacterial load.

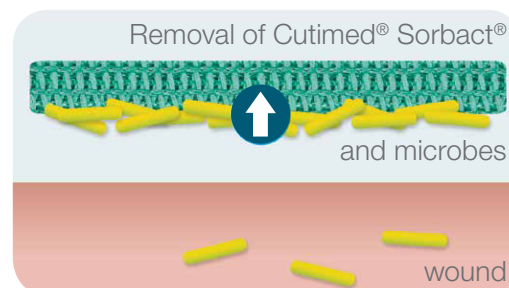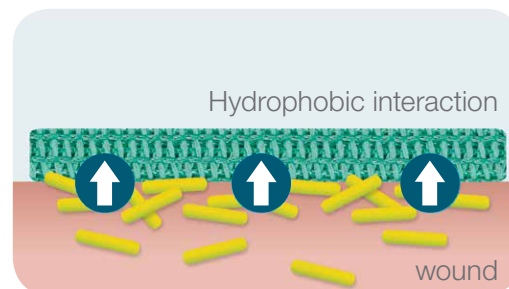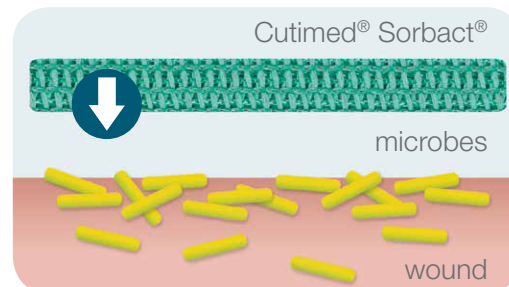

Wound Depth	Wound Phase	Exudate
<input type="checkbox"/> Superficial	<input type="checkbox"/> Necrotic	<input type="checkbox"/> Low
<input type="checkbox"/> Deep	<input type="checkbox"/> Infected	<input type="checkbox"/> Low to moderate
<input type="checkbox"/> Superficial + deep	<input type="checkbox"/> Sloughy	<input type="checkbox"/> Moderate to high
	<input type="checkbox"/> Granulation	<input type="checkbox"/> High to Excessive
	<input type="checkbox"/> Epithelialization	

Description		Dimensions	Code	Box Contents	Selling Unit
Cutimed® Sorbact® pad With absorbent core, ideal under TCC 		7 cm x 9 cm	7216101	5 dressings	Box
		7 cm x 9 cm	7216100	40 dressings	Box
		10 cm x 10 cm	7216201	5 dressings	Box
		10 cm x 10 cm	7216200	40 dressings	Box
		10 cm x 20 cm	7216300	20 dressings	Box
Cutimed® Sorbact® swab Most versatile; for any exudate level, deep or superficial wounds 		4 cm x 6 cm (11 x 16 unfolded)	7216401	5 dressings	Box
		4 cm x 6 cm (11 x 16 unfolded)	7216400	40 dressings	Box
		7 cm x 9 cm (17 x 28 unfolded)	7216501	5 dressings	Box
		7 cm x 9 cm (17 x 28 unfolded)	7216500	40 dressings	Box
Cutimed® Sorbact® ribbon For use in cavity wounds of all kinds 		2 cm x 50 cm	7216600	20 dressings	Box
		5 cm x 200 cm	7216700	10 dressings	Box
Cutimed® Sorbact® round swab For deep wound dressing or superficial swabing 		3.5 cm (round)	7216800	70 dressings (14 x 5)	Box
Cutimed® Sorbact® gel Impregnated with hydrogel; for dry wounds 		7.5 cm x 7.5 cm	7261100	10 dressings	Box
		7.5 cm x 15 cm	7261101	10 dressings	Box
Cutimed® Sorbact® Hydroactive Absorbent hydrogel core 		7 cm x 8.5 cm	7264603	10 dressings	Box
		14 cm x 14 cm	7264604	10 dressings	Box
		19 cm x 19 cm	7264605	10 dressings	Box
		14 cm x 24 cm	7264607	10 dressings	Box
		24 cm x 24 cm	7264606	5 dressings	Box
Cutimed® Sorbact® Hydroactive B Absorbent hydrogel core and skin-friendly adhesive border 		7 cm x 8.5 cm	7993300	10 dressings	Box
		14 cm x 14 cm	7993302	10 dressings	Box
		19 cm x 19 cm	7993303	10 dressings	Box
		14 cm x 24 cm	7993304	10 dressings	Box
Cutimed® Siltec Sorbact® Super-absorbent foam dressing with silicone adhesive border 		7.5 cm x 7.5 cm	7992900	10 dressings	Box
		12.5 cm x 12.5 cm	7992901	10 dressings	Box
		15 cm x 15 cm	7992902	10 dressings	Box
		17.5 cm x 17.5 cm	7992903	5 dressings	Box
		22.5 cm x 22.5 cm	7992904	5 dressings	Box
		Sacrum 17.5 cm x 17.5 cm	7992905	5 dressings	Box
		Sacrum 23 cm x 23 cm	7992906	5 dressings	Box
Cutimed® Sorbion® Sorbact® Super-absorbent core with gelling effect to retain very high levels of exudate 		10 cm x 10 cm	7269800	10 dressings	Box
		10 cm x 20 cm	7269801	10 dressings	Box
		20 cm x 20 cm	7269802	10 dressings	Box
		20 cm x 30 cm	7269803	10 dressings	Box

CUTIMED® SILTEC

FOAM DRESSING WITH SILICONE WOUND CONTACT LAYER



Cutimed® Siltec consistently provides reliable results, regardless of exudate level or viscosity. With its silicone wound contact layer and super-absorbent particles, it offers outstanding fluid management capabilities, but also offers patients key benefits for increased quality of life and compliance.

INDICATIONS

- For all exuding wounds such as venous, diabetic or pressure ulcers

FEATURES AND BENEFITS

- Reduced risk of maceration
- Atraumatic, pain-free dressing changes
- Visible saturation
- Strikethrough and leakage prevention
- Comfortable wear
- Icons on box to choose the right product
- Sterile
- Latex Free

Description	Dimensions	Code	Box Contents	Selling Unit
Cutimed® Siltec With very gentle <i>FeatherTack</i> adherence 	 5 cm x 6 cm	7263200	10 dressings	Box
	10 cm x 10 cm	7263201	10 dressings	Box
	10 cm x 20 cm	7263202	10 dressings	Box
	15 cm x 15 cm	7263203	10 dressings	Box
	20 cm x 20 cm	7263204	5 dressings	Box
	Heel 16 cm x 24 cm	7264800	5 dressings	Box
Cutimed® Siltec^{PLUS} With <i>SoftTack</i> adherence  	5 cm x 6 cm	7268300	10 dressings	Box
	10 cm x 10 cm	7268301	10 dressings	Box
	10 cm x 20 cm	7268302	10 dressings	Box
	15 cm x 15 cm	7268303	10 dressings	Box
	20 cm x 20 cm	7268304	5 dressings	Box
	Heel 16 cm x 24 cm	7268305	5 dressings	Box
Cutimed® Siltec B (Border) With silicone adhesive border 	 7.5 cm x 7.5 cm	7263100	10 dressings	Box
	10 cm x 10 cm	7263105	10 dressings	Box
	10 cm x 22.5 cm	7263106	10 dressings	Box
	12.5 cm x 12.5 cm	7263101	10 dressings	Box
	15 cm x 15 cm	7263102	10 dressings	Box
	17.5 cm x 17.5 cm	7263103	5 dressings	Box
	22.5 cm x 22.5 cm	7263104	5 dressings	Box
Cutimed® Siltec Sacrum For sacral wounds 	 17.5 cm x 17.5 cm	7264700	5 dressings	Box
	23 cm x 23 cm	7264701	5 dressings	Box
Cutimed® Siltec L (Light) For low exuding wounds 	5 cm x 6 cm	7263000	10 dressings	Box
	10 cm x 10 cm	7263001	10 dressings	Box
	15 cm x 15 cm	7263002	10 dressings	Box
Cutimed® Siltec Sorbact® Combined with Sorbact® antimicrobial technology; with skin-friendly silicone adhesive border 	7.5 cm x 7.5 cm	7992900	10 dressings	Box
	12.5 cm x 12.5 cm	7992901	10 dressings	Box
	15 cm x 15 cm	7992902	10 dressings	Box
	17.5 cm x 17.5 cm	7992903	5 dressings	Box
	22.5 cm x 22.5 cm	7992904	5 dressings	Box
	Sacrum 17.5 cm x 17.5 cm	7992905	5 dressings	Box
	Sacrum 23 cm x 23 cm	7992906	5 dressings	Box

* See legend on page 18

CUTIMED® SORBION®

SUPER-ABSORBENT WOUND DRESSING

Cutimed® Sorbion® dressings have the excellent ability of quickly absorbing large amounts of exudate and locking it within their inner core, even under compression, therefore reducing the frequency of dressing changes. Cutimed® Sorbion® dressings also have the potential to remove chronicity-causing factors and thereby support wound healing; they provide gentle debridement, reduce the concentration and activity of MMPs and immobilize micro-organisms. The unique fiber matrix offers stability even when highly saturated, thus preventing the formation of lumps and pressure points and ensuring maximum wound contact.

INDICATIONS

- For all moderately to highly exuding wounds such as venous, diabetic or pressure ulcers

FEATURES AND BENEFITS

- High absorption capacity, long dressing wear time
- Fewer dressing changes (up to 4 days wear time)
- Securely locks in exudate, even under compression
- Removes wound debris and locks in wound germs
- Maintains a moist wound environment by gel formation
- No binders or adhesives; well tolerated
- No risk of leakage, high tear strength
- Able to apply on both sides
- Less odor formation
- Sterile
- Latex Free



Description		Dimensions	Code	Box Contents	Selling Unit
Cutimed® Sorbion® Sachet S The classic form suited for many wound types and areas of application 		7.5 cm x 7.5 cm	7323200	10 dressings	Box
		10 cm x 10 cm	7323206	10 dressings	Box
		5 cm x 12 cm	7323203	10 dressings	Box
		10 cm x 20 cm	7323209	10 dressings	Box
		15 cm x 15 cm	7323212	10 dressings	Box
		20 cm x 20 cm	7323215	10 dressings	Box
		20 cm x 30 cm	7323218	10 dressings	Box
Cutimed® Sorbion® Sachet Multi Star Flexible form for applications on toes, heel, foot, sacrum, chest or wound cavities 		8 cm	7323701	10 dressings	Box
		14 cm	7323703	10 dressings	Box
Cutimed® Sorbion® Sachet XL For large difficult-to-dress wounds on lower leg, sacral region, thorax, thigh, lower part of abdomen 		25 cm x 45 cm	7324000	10 dressings	Box
Cutimed® Sorbion® Sachet Border With hypoallergene adhesive border 		10 cm x 10 cm	7323600	10 dressings	Box
		15 cm x 15 cm	7323602	10 dressings	Box
		15 cm x 25 cm	7323604	10 dressings	Box
Cutimed® Sorbion® Sana With particularly atraumatic surface for pain-sensitive patients 		8.5 cm x 8.5 cm	7323300	10 dressings	Box
		12 cm x 12 cm	7323303	10 dressings	Box
		12 cm x 22 cm	7323306	10 dressings	Box
		22 cm x 22 cm	7323309	10 dressings	Box
Cutimed® Sorbion® Sana Multi Star With particularly atraumatic surface for pain-sensitive patients 		11 cm	7323801	10 dressings	Box
		17 cm	7323803	10 dressings	Box

* See legend on page 18

CUTIMED® CAVITY

ABSORBENT FOAM DRESSING FOR DEEP WOUNDS



Cutimed® Cavity is highly absorbent and conformable, making it the ideal choice for deep exuding wounds.

INDICATIONS

- For all deep exuding wounds such as venous, diabetic or pressure ulcers
- Wounds after excision of fistulas and abscesses

FEATURES AND BENEFITS

- Soft and conformable
- Atraumatic
- Non-adhesive
- Easy to apply and remove
- Highly absorbent
- Icons on box to choose the right product*
- Sterile
- Latex Free



Dimensions	Code	Box Contents	Selling Unit
5 cm x 6 cm	7262100	10 dressings	Box
10 cm x 10 cm	7262101	10 dressings	Box
15 cm x 15 cm	7262103	5 dressings	Box

* See legend on page 18

CUTIMED® TOTAL CONTACT CAST KIT

PRESSURE OFFLOADING CAST



The Cutimed® Total Contact Cast Kit has been developed in partnership with healthcare professionals dealing with diabetic foot ulcers every day. The result is a recipe that offers the most effective TCC through an easy standardized technique. It combines specifically chosen and proven materials to provide an intimate close fit and proper pressure redistribution for an optimized healing environment.

INDICATIONS

- Diabetic foot ulcers
- Charcot joint injuries
- Post-op immobilization following deformity reconstructions

FEATURES AND BENEFITS

- Increased patient compliance
- Optimized healing environment
- Standardized technique
- Convenient all-in-one kit
- Cost-effective
- Latex Free

Code	Box Contents	Selling Unit
7349900	All Material for 1 Cast	Box

LEUKOSAN® STRIP

ADHESIVE CLOSURE STRIPS



Leukosan® Strip's prevalent benefits include atraumatic, painless application, low infection risk and good cosmetic results.

INDICATIONS

- Facial lacerations or incisions
- Thin, fragile or sensitive skin
- Wounds closed with intracutaneous sutures
- Replacement of staples and sutures

FEATURES AND BENEFITS

- Atraumatic, pain-free application
- Tension-free adaptation of wound edges
- Good cosmetic results
- Low infection risk
- Excellent tolerance
- Breathable
- Sterile
- Latex Free

Dimensions	Code	Box Contents	Case Contents	Selling Unit
6 mm x 38 mm	7262900	50 x 6 strips	4 boxes	Case
3 mm x 75 mm	7262901	50 x 5 strips	4 boxes	Case
6 mm x 75 mm	7262902	50 x 3 strips	4 boxes	Case
6 mm x 100 mm	7262903	50 x 10 strips	4 boxes	Case
12 mm x 100 mm	7262904	50 x 6 strips	4 boxes	Case
25 mm x 100 mm	7262905	25 x 4 strips	4 boxes	Case

CUTIMED® GEL

TRANSPARENT AMORPHOUS HYDROGEL



Cutimed® Gel is a clear, amorphous hydrogel with pure formulation to reduce the risk of allergic reaction, and its superior moisture donation allows for less frequent dressing changes and effective debridement.

INDICATIONS

- To support autolytic debridement in necrotic wounds
- To provide moisture to dry wounds

FEATURES AND BENEFITS

- No preservatives
- Icons on box to choose the right product*
- Sterile
- Latex Free



Dimensions	Code	Box Contents	Selling Unit
Tube 8 g	7261000	10 tubes	Box
Tube 15 g	7261001	10 tubes	Box
Tube 25 g	7261002	10 tubes	Box

* See legend on page 18

COMPRILAN®

SHORT STRETCH COMPRESSION BANDAGE



Unlike long stretch bandages, Comprilan® short stretch provides high working pressure and low resting pressure, making it better tolerated by patients. It enables for cost-effective treatment with easy reapplication.

INDICATIONS

- Venous and mixed leg ulcers
- Edema
- Lymphedema
- Thrombophlebitis
- Post sclerotherapy

FEATURES AND BENEFITS

- Better tolerated by patients
- Safe for mixed ulcers ($0.6 < \text{ABPI}$)
- Cost-effective when frequent dressing changes required (washable and reusable up to 10 times)
- Easy to apply and safe to use, even with minimal training
- As effective as multi-layer bandaging
- Latex Free (except clips)

Dimensions	Code	Box Contents	Selling Unit
4 cm x 5 m	7718700	1 roll	Box
6 cm x 5 m	0102600	1 roll	Box
8 cm x 5 m	0102700	1 roll	Box
10 cm x 5 m	0102800	1 roll	Box
12 cm x 5 m	0102900	1 roll	Box
10 cm x 10 m	7718800	1 roll	Box
12 cm x 10 m	7718900	1 roll	Box



JOBST® COMPRI2

2-LAYER BANDAGE SYSTEM



JOBST® Compriz, the new 2-layer short stretch compression system, addresses the two key elements required to ensure successful compression therapy: correct compression level easily applied and continuous patient adherence.

INDICATIONS

- Venous and mixed leg ulcers
- Lymphedema

FEATURES AND BENEFITS

- Provides sustained graduated compression of 40 mmHg (JOBST® Compriz2) or 20-30 mmHg (JOBST® Compriz2 Lite) for up to 7 days
- Easy to apply with 2-way stretch padding and compression indicators
- Corrective action / re-application possible
- Comfortable skin-friendly padding
- Thin system
- Online support (application videos, micro-site for patients and professionals)
- Latex Free

Description	Code	Ankle Size	Box Contents	Selling Unit
JOBST® Compriz2	7627100	18-25 cm	1 kit	Box
	7627101	25-32 cm	1 kit	Box
JOBST® Compriz2 lite	7627102	18-25 cm	1 kit	Box
	7627103	25-32 cm	1 kit	Box

JOBST® COMPRIFORE LF

MULTI-LAYER COMPRESSION BANDAGE



JOBST® Comprifore LF is a multi-layer compression system for the treatment of leg ulcers. With the included Cutimed® Sorbact® WCL (wound contact layer), the set also allows for management of bacterial burden.

INDICATIONS

- Venous and mixed leg ulcers

FEATURES AND BENEFITS

- Provides sustained graduated compression of 40 mmHg (JOBST® Comprifore LF) or 20-30 mmHg (JOBST® Comprifore Lite LF) for up to 7 days
- Contains Cutimed® Sorbact® WCL, with bacteria-binding action and low adherence for atraumatic and pain-free dressing changes
- Latex Free

Description	Code	Box Contents	Selling Unit
JOBST® Comprifore LF (4 Layers)	7266101	1 kit	Box
JOBST® Comprifore Lite LF (3 Layers)	7266103	1 kit	Box

TENSOPRESS®

LONG STRETCH COMPRESSION BANDAGE



Tensopress® delivers measured, sustained compression for the treatment of venous leg ulcers.

INDICATIONS

- Venous leg ulcers
- Venous insufficiency
- Deep vein thrombosis
- Compression of stumps

FEATURES AND BENEFITS

- Soft edges prevent cutting in
- Washable
- Provides sustained, even and consistent 40 mmHg compression for up to 8 days
- Central green line which guides and simplifies precise and even application
- Long stretch

Dimensions	Code	Box Contents	Case Contents	Selling Unit
7.5 cm x 3 m	7172301	1 roll	10 boxes	Case
10 cm x 3 m	7172300	1 roll	10 boxes	Case

GELOCAST®

ZINC PASTE COMPRESSION BANDAGE (UNNA BOOT)



Gelocast® the original Unna boot is a gauze bandage impregnated with zinc oxide and calamine. Used as an aid in the management of venous leg ulcers, Gelocast® provides high working pressure and lower resting pressure to reduce edema.

INDICATIONS

- Venous and mixed leg ulcers
- Phlebitis
- Strains and sprains
- Dislocations
- As an after-cast

FEATURES AND BENEFITS

- Dries quickly providing a semi-rigid cast
- 100% cotton base does not unravel reducing wastage
- Provides high working pressure and lower resting pressure
- Uniformly impregnated, easy to apply and remove
- Latex Free

Dimensions	Code	Box Contents	Selling Unit
7.5 cm x 9 m	105200	1 roll	Box
10 cm x 9 m	105300	1 roll	Box

JOBST® ULCERCARE™

READY-TO-WEAR COMPRESSION SYSTEM



JOBST® UlcerCARE™ enables effective management of a venous leg ulcer while offering ease of use at the same time – the perfect combination for providing patients with a higher quality of life.

INDICATIONS

- Venous leg ulcers

FEATURES AND BENEFITS

- Provides graduated compression of 40 mmHg
- High stiffness factor for effective compression
- Easy handling by patients or caregivers
- Contains silk in the liner for ease of donning
- Supports patients' independence
- Compatible with common footwear and clothing to enable a normal lifestyle
- Latex Free

Style	Colour	S	M	L	XL	2XL	3XL	4XL	Box Contents	Selling Unit
No Zipper	Beige	7363021	7363022	7363023	7363024	7363025	7363026	7363027	1 stocking & 2 liners	Box
No Zipper	Black	7363041	7363042	7363043	7363044	7363045	7363046	7363047	1 stocking & 2 liners	Box
Left Zipper	Beige	7363121	7363122	7363123	7363124	7363125	7363126	7363127	1 stocking & 2 liners	Box
Left Zipper	Black	7363141	7363142	7363143	7363144	7363145	7363146	7363147	1 stocking & 2 liners	Box
Right Zipper	Beige	7363161	7363162	7363163	7363164	7363165	7363166	7363167	1 stocking & 2 liners	Box
Right Zipper	Black	7363171	7363172	7363173	7363174	7363175	7363176	7363177	1 stocking & 2 liners	Box
Replacement Liners										
	White	7363221	7363222	7363223	7363224	7363225	7363226	7363227	3 liners	Box

SIZING CHART

Measurements (cm)	S	M	L	XL	2XL	3XL	4XL
Ankle	18-21	21-24	24-27	27-30	30-33	33-36	36-39
Calf	29-35	34-41	40-47	46-54	48-56	50-58	52-60

JOBST® FARROWWRAP®

SHORT STRETCH COMPRESSION WRAP



JOBST® FarrowWrap® provides steady and comfortable resting compression, but higher working compression when the patient is active. It helps prevent limb swelling during the day. Since it is applied at full-stretch, the more the limb tries to swell, the more the garment works to prevent swelling.

INDICATIONS

- Chronic venous insufficiency
- Mild to moderate edema
- Lymphedema

FEATURES AND BENEFITS

- Easy to don for patients with limited mobility
- Easy to use adjustable bands
- Supports patients' independence
- Available for lower and upper extremity in ready-to-wear, trim-to-fit and custom
- Latex Free

SIZING CHART

Measurement (cm)	XS	S	M	L	XL
Thighpiece					
Groin circumference	60-70	67-75	70-80	75-90	85-100
Above knee circumference	47-52	53-59	60-66	67-73	74-80
E-G Short			16-20		
E-G Regular			21-26		
E-G Long			27-31		
Legpiece (except JOBST® FarrowWrap® 4000)					
Mid-calf circumference	36-43	42-50	48-58	53-63	58-68
Ankle circumference	21-25	25-30	30-36	36-42	42-50
A-D Regular*	33-37	35-39	37-41	39-43	39-43
A-D Long*	38-41	40-43	42-45	44-47	44-47
B-D Regular**	28	30	31	33	33
B-D Long**	33	35	36	38	38
Legpiece JOBST® FarrowWrap® 4000					
Mid-calf circumference	25-38	32-48	37-56	43-67	-
Ankle circumference	17-23	19-28	25-37	29-43	-
B-D Regular**			32		
B-D Long**			36		
Footpiece					
Midfoot circumference	22-24	25-27	28-30	31-34	35-40
Regular	16-17	17-18	19-20	20-21	22-23
Long	18-19	19-20	21-22	22-23	24-25
Liners***					
TG® Soft (calf circumference)	-	<40	40-70	70-120	-
Flesh Tone Silver Liners (calf circumference)	-	19-48	-	48-70	-

*From floor to 2 fingers below knee crease

**B-D: Actual height of the garment

***Supplied in pairs

ACCESSORIES

Description	Colour	S	M	L	XL
JOBST® FarrowHybrid® Foot Compression					
ADI Regular 20-30 mmHg	Tan	7666400	7666401	7666402	7666403
ADI Wide 20-30 mmHg	Tan	7666412	7666413	7666414	7666415
ADI Wide 20-30 mmHg	Black	7666416	7666417	7666418	7666419
ADII Regular 30-40 mmHg	Tan	7666404	7666405	7666406	7666407
ADII Wide 30-40 mmHg	Tan	7666408	7666409	7666410	7666411
TG® Soft Liners					
Knee High Open Toe	Tan	7666900	7666902	7666904	-
Thigh High Open Toe	Tan	7667000	7667002	7667004	-
Flesh Tone Silver Liners					
Knee High	Tan	7666800	-	7666802	-

Description	Colour	XS	S	M	L	XL
Thighpiece						
LITE						
Short	Tan	7665505	7665506	7665507	7665508	7665509
Regular	Tan	7665500	7665501	7665502	7665503	7665504
Long	Tan	7665510	7665511	7665512	7665513	7665514
CLASSIC						
Short	Tan	7665806	7665807	7665808	7665809	7665810
Regular	Tan	7665800	7665801	7665802	7665803	7665804
Long	Tan	7665812	7665813	7665814	7665815	7665816
STRONG						
Short	Tan	7664905	7664906	7664907	7664908	7664909
Regular	Tan	7664900	7664901	7664902	7664903	7664904
Long	Tan	7664910	7664911	7664912	7664913	7664914
Legpiece						
LITE						
Regular	Tan	7665400	7665401	7665402	7665403	7665404
Long	Tan	7665405	7665406	7665407	7665408	7665409
CLASSIC						
Regular	Tan	7665700	7665701	7665702	7665703	7665704
Long	Tan	7665705	7665706	7665707	7665708	7665709
STRONG						
Regular	Tan	7664800	7664801	7664802	7664803	7664804
Long	Tan	7664805	7664806	7664807	7664808	7664809
BASIC*						
Regular (includes Velcro®)	Tan	7666000	7666001	7666002	7666003	7666004
Long (includes Velcro®)	Tan	7666005	7666006	7666007	7666008	7666009
Extra Velcro®	7666010					
4000						
Regular	Tan	7666100	7666101	7666102	7666103	-
Long	Tan	7666104	7666105	7666106	7666107	-
Regular	Black	7666108	7666109	7666110	7666111	-
Long	Black	7666112	7666113	7666114	7666115	-
Footpiece						
LITE						
Regular	Tan	7665205	7665206	7665207	7665208	7665209
Long	Tan	7665200	7665201	7665202	7665203	7665204
CLASSIC						
Regular	Tan	7665605	7665606	7665607	7665608	7665609
Long	Tan	7665600	7665601	7665602	7665603	7665604
STRONG						
Regular	Tan	7664604	7664605	7664606	7664607	7664608
Long	Tan	7664609	7664600	7664601	7664602	7664603
BASIC*						
Regular	Tan	7665905	7665906	7665907	7665908	7665909
Long	Tan	7665900	7665901	7665902	7665903	7665904

*Free pair of JOBST® FarrowHybrid® ADI Wide Foot Compression Liners included

COMPRIFOAM®

FOAM BANDAGE



Comprifoam® open cell foam bandage is designed to evenly distribute the pressure from compression bandages over extremities, thus preventing constrictions and ensuring an effective compression therapy.

INDICATIONS

- Padding under compression bandages

FEATURES AND BENEFITS

- Breathable
- Non-sliding
- Tear resistant
- Easy for patients to self bandage
- Anatomically conforming allowing freedom of movement
- Resistant to body oils, lotions and perspiration
- Reusable
- May be steam sterilized
- Latex Free

Dimensions	Code	Box Contents	Selling Unit
10 cm x 2.5 m x 0.4 cm	7529400	1 roll	Box
12 cm x 2.5 m x 0.4 cm	7529500	1 roll	Box

ARTIFLEX®

NON-WOVEN PADDING BANDAGE



Artiflex® is a multipurpose padding bandage that ensures even distribution of pressure in phlebological or lymphedema applications. Also useful as under padding of plaster or synthetic casts.

INDICATIONS

- Padding under compression bandages for sensitive skin
- Under padding of plaster and synthetic casts

FEATURES AND BENEFITS

- Non-absorbent
- Applies easily without creasing
- Provides added comfort in wrapping for lymphedema
- Breathable
- Autoclavable
- Latex Free

Dimensions	Code	Box Contents	Selling Unit
10 cm x 3 m	09046	30 rolls	Box
15 cm x 3 m	09047	20 rolls	Box

FOAM RUBBER PAD



Dimensions	Code	Box Contents	Selling Unit
9 cm x 1 cm	7849600	Unit	Box

ANTI-EM/GP™

ANTI-EMBOLISM STOCKINGS WITH GRADUATED COMPRESSION



JOBST® Anti-Em/GP™ anti-embolism stockings are specifically designed to combat thromboembolism. They apply graduated compression, greatest at the ankle and gradually decreasing up the leg without tourniqueting or impeding the venous return.


INDICATIONS

- Prevention of thromboembolism

FEATURES AND BENEFITS

- Inspection window on top of the foot so the patients' skin doesn't touch the floor
- Toe colour coordinated with sizing chart
- Latex Free

SIZING CHARTS

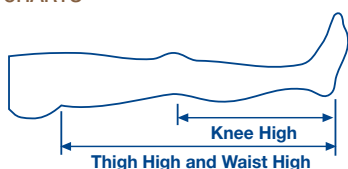

Measurements (cm)	S	M	L	XL	2XL
Calf	23-30	30-38	38-44	44-51	46-53
Thigh	41-74	41-74	41-74	61-79	81-90
Knee high – Regular			36-43		
Knee high – Long			43-51		
Thigh high – Short			64-74		
Thigh/waist high – Regular			74-84		
Thigh/waist high – Long			84-97		

Size	Code	Box Contents	Selling Unit
Knee High			
Small Regular	111402	12 pairs	Box
Small Long	111403	12 pairs	Box
Medium Regular	111406	12 pairs	Box
Medium Long	111407	12 pairs	Box
Large Regular	111410	12 pairs	Box
Large Long	111411	12 pairs	Box
XLarge Regular	111414	12 pairs	Box
XLarge Long	111415	12 pairs	Box
Thigh High			
Small Short	111450	6 pairs	Box
Small Regular	111451	6 pairs	Box
Small Long	111452	6 pairs	Box
Medium Short	111454	6 pairs	Box
Medium Regular	111455	6 pairs	Box
Medium Long	111456	6 pairs	Box
Large Short	111458	6 pairs	Box
Large Regular	111459	6 pairs	Box
Large Long	111460	6 pairs	Box
XLarge Regular	111462	6 pairs	Box
XLarge Long	111463	6 pairs	Box
2XLarge Regular	111464	6 pairs	Box
2XLarge Long	111465	6 pairs	Box
Waist High			
Small Regular	111625	6 pairs	Box
Small Long	111626	6 pairs	Box
Medium Regular	111627	6 pairs	Box
Medium Long	111628	6 pairs	Box
Large Regular	111629	6 pairs	Box
Large Long	111630	6 pairs	Box

DELTA-LITE® PLUS

FIBERGLASS CAST TAPE



Delta-Lite® Plus provides extra rigid casts with a smooth surface and edges. The advanced water activated resin system provides tack free properties allowing for ease of application.

FEATURES AND BENEFITS

- Smooth surface
- Rigid support
- Good conformability
- Excellent lamination
- Outstanding molding
- Latex Free

Dimensions	Colour	Code	Box Contents	Selling Unit
2.5 cm x 1.8 m	White	7345800	10 rolls	Box
5 cm x 3.6 m	White	7345801	10 rolls	Box
7.5 cm x 3.6 m	White	7345802	10 rolls	Box
10 cm x 3.6 m	White	7345803	10 rolls	Box
12.5 cm x 3.6 m	White	7345804	10 rolls	Box
5 cm x 3.6 m	Dark Blue	7345820	10 rolls	Box
7.5 cm x 3.6 m	Dark Blue	7345821	10 rolls	Box
10 cm x 3.6 m	Dark Blue	7345822	10 rolls	Box
12.5 cm x 3.6 m	Dark Blue	7345823	10 rolls	Box
5 cm x 3.6 m	Green	7345825	10 rolls	Box
7.5 cm x 3.6 m	Green	7345826	10 rolls	Box
10 cm x 3.6 m	Green	7345827	10 rolls	Box
12.5 cm x 3.6 m	Green	7345828	10 rolls	Box
5 cm x 3.6 m	Red	7345830	10 rolls	Box
7.5 cm x 3.6 m	Red	7345831	10 rolls	Box
10 cm x 3.6 m	Red	7345832	10 rolls	Box
5 cm x 3.6 m	Light Blue	7345835	10 rolls	Box
7.5 cm x 3.6 m	Light Blue	7345836	10 rolls	Box
10 cm x 3.6 m	Light Blue	7345837	10 rolls	Box
12.5 cm x 3.6 m	Light Blue	7345838	10 rolls	Box
5 cm x 3.6 m	Black	7345845	10 rolls	Box
7.5 cm x 3.6 m	Black	7345846	10 rolls	Box
10 cm x 3.6 m	Black	7345847	10 rolls	Box
5 cm x 3.6 m	Mixed Colours*	7345850	10 rolls	Box
7.5 cm x 3.6 m	Mixed Colours*	7345851	10 rolls	Box
10 cm x 3.6 m	Mixed Colours*	7345852	10 rolls	Box
5 cm x 3.6 m	Pink	7345855	10 rolls	Box
7.5 cm x 3.6 m	Pink	7345856	10 rolls	Box
10 cm x 3.6 m	Pink	7345857	10 rolls	Box
5 cm x 3.6 m	Purple	7345860	10 rolls	Box
7.5 cm x 3.6 m	Purple	7345861	10 rolls	Box
10 cm x 3.6 m	Purple	7345862	10 rolls	Box
5 cm x 3.6 m	Yellow	7345865	10 rolls	Box
7.5 cm x 3.6 m	Yellow	7345866	10 rolls	Box
10 cm x 3.6 m	Yellow	7345867	10 rolls	Box
Reinforcing Slabs				
5 cm x 20 cm	White	7348500	10 splints	Box
7.5 cm x 20 cm	White	7348501	10 splints	Box
7.5 cm x 38 cm	White	7227200	10 splints	Box
10 cm x 38 cm	White	7348503	10 splints	Box
7.5 cm x 76 cm	White	7348502	10 splints	Box
10 cm x 90 cm	White	7227201	10 splints	Box
12.5 cm x 90 cm	White	7348504	10 splints	Box

* Mixed Colours: Dark Blue, Red, Light Blue, Purple and Pink

DELTA-CAST® CONFORMABLE

POLYESTER CAST TAPE



Delta-Cast® Conformable is an easy-to-use polyester cast tape that produces versatile, lightweight and strong casts with unsurpassed conformability. Ideal for applications where patient fit and comfort are critical.

FEATURES AND BENEFITS

- Unsurpassed conformability
- Smooth finish and soft edges
- Rigid support
- 100% radiolucent
- Excellent lamination
- Univalve
- Latex Free

Dimensions	Colour	Code	Box Contents	Selling Unit
2.5 cm x 1.8 m	White	7228000	10 rolls	Box
5 cm x 3.6 m	White	7228001	10 rolls	Box
7.5 cm x 3.6 m	White	7228002	10 rolls	Box
10 cm x 3.6 m	White	7228003	10 rolls	Box
12.5 cm x 3.6 m	White	7228004	10 rolls	Box
5 cm x 3.6 m	Dark Blue	7228200	10 rolls	Box
7.5 cm x 3.6 m	Dark Blue	7228201	10 rolls	Box
10 cm x 3.6 m	Dark Blue	7228202	10 rolls	Box
5 cm x 3.6 m	Green	7228700	10 rolls	Box
7.5 cm x 3.6 m	Green	7228701	10 rolls	Box
10 cm x 3.6 m	Green	7228702	10 rolls	Box
5 cm x 3.6 m	Red	7228500	10 rolls	Box
7.5 cm x 3.6 m	Red	7228501	10 rolls	Box
10 cm x 3.6 m	Red	7228502	10 rolls	Box
5 cm x 3.6 m	Black	7228800	10 rolls	Box
7.5 cm x 3.6 m	Black	7228801	10 rolls	Box
10 cm x 3.6 m	Black	7228802	10 rolls	Box
5 cm x 3.6 m	Mixed Colours*	7228910	10 rolls	Box
7.5 cm x 3.6 m	Mixed Colours*	7228911	10 rolls	Box
10 cm x 3.6 m	Mixed Colours*	7228912	10 rolls	Box
5 cm x 3.6 m	Pink	7228100	10 rolls	Box
7.5 cm x 3.6 m	Pink	7228101	10 rolls	Box
10 cm x 3.6 m	Pink	7228102	10 rolls	Box
5 cm x 3.6 m	Purple	7228600	10 rolls	Box
7.5 cm x 3.6 m	Purple	7228601	10 rolls	Box
10 cm x 3.6 m	Purple	7228602	10 rolls	Box
5 cm x 3.6 m	Yellow	7228300	10 rolls	Box
7.5 cm x 3.6 m	Yellow	7228301	10 rolls	Box
10 cm x 3.6 m	Yellow	7228302	10 rolls	Box
5 cm x 3.6 m	Orange	7228400	10 rolls	Box
7.5 cm x 3.6 m	Orange	7228401	10 rolls	Box
10 cm x 3.6 m	Orange	7228402	10 rolls	Box

* Mixed Colours: Dark Blue, Red, Black, Purple and Pink

DELTA-CAST® ELITE

POLYESTER CAST TAPE



Delta-Cast® Elite is a polyester cast tape that offers strength and durability. It features extensible yarns that provide excellent molding capabilities and a smooth cast surface for patient comfort.

FEATURES AND BENEFITS

- Smooth finish and soft edges
- Rigid support
- 100% radiolucent
- Excellent conformability
- Excellent lamination
- Univalve
- Latex Free

Dimensions	Colour	Code	Box Contents	Selling Unit
2.5 cm x 1.8 m	White	7227904	10 rolls	Box
5 cm x 3.6 m	White	7227900	10 rolls	Box
7.5 cm x 3.6 m	White	7227901	10 rolls	Box
10 cm x 3.6 m	White	7227902	10 rolls	Box
12.5 cm x 3.6 m	White	7227903	10 rolls	Box
5 cm x 3.6 m	Dark Blue	7227600	10 rolls	Box
7.5 cm x 3.6 m	Dark Blue	7227601	10 rolls	Box
10 cm x 3.6 m	Dark Blue	7227602	10 rolls	Box
12.5 cm x 3.6 m	Dark Blue	7227603	10 rolls	Box
5 cm x 3.6 m	Green	7345400	10 rolls	Box
7.5 cm x 3.6 m	Green	7345401	10 rolls	Box
10 cm x 3.6 m	Green	7345402	10 rolls	Box
5 cm x 3.6 m	Red	7227500	10 rolls	Box
7.5 cm x 3.6 m	Red	7227501	10 rolls	Box
10 cm x 3.6 m	Red	7227502	10 rolls	Box
5 cm x 3.6 m	Light Blue	7345300	10 rolls	Box
7.5 cm x 3.6 m	Light Blue	7345301	10 rolls	Box
10 cm x 3.6 m	Light Blue	7345302	10 rolls	Box
5 cm x 3.6 m	Black	7227800	10 rolls	Box
7.5 cm x 3.6 m	Black	7227801	10 rolls	Box
10 cm x 3.6 m	Black	7227802	10 rolls	Box
5 cm x 3.6 m	Mixed Colours*	7345700	10 rolls	Box
7.5 cm x 3.6 m	Mixed Colours*	7345701	10 rolls	Box
10 cm x 3.6 m	Mixed Colours*	7345702	10 rolls	Box
5 cm x 3.6 m	Pink	7345500	10 rolls	Box
7.5 cm x 3.6 m	Pink	7345501	10 rolls	Box
10 cm x 3.6 m	Pink	7345502	10 rolls	Box
5 cm x 3.6 m	Purple	7227700	10 rolls	Box
7.5 cm x 3.6 m	Purple	7227701	10 rolls	Box
10 cm x 3.6 m	Purple	7227702	10 rolls	Box
5 cm x 3.6 m	Orange	7345600	10 rolls	Box
7.5 cm x 3.6 m	Orange	7345601	10 rolls	Box
10 cm x 3.6 m	Orange	7345602	10 rolls	Box
Reinforcing Slabs 3 Ply				
5 cm x 20 cm	White	7454507	10 splints	Box
7.5 cm x 20 cm	White	7454508	10 splints	Box
7.5 cm x 38 cm	White	7454509	10 splints	Box
Reinforcing Slabs 5 Ply				
10 cm x 38 cm	White	7454503	10 splints	Box
10 cm x 90 cm	White	7454504	10 splints	Box
12.5 cm x 90 cm	White	7454505	10 splints	Box

* Mixed Colours: Dark Blue, Red, Light Blue, Purple and Pink

DELTA-CAST® PRINTS

POLYESTER PRINTED CAST TAPE



Add personality to casts with colourful printed casting material. In addition to helping relieve patient anxiety and making casts more enjoyable to wear, the finished cast is strong, lightweight and durable. Patterns are colourful and will not interfere with radiographic viewing.

FEATURES AND BENEFITS

- Same advantages as Delta-Cast® Elite
- Unique style
- Long lasting, vibrant prints
- Latex Free

Dimensions	Colour	Code	Box Contents	Selling Unit
5 cm x 3.6 m	Pastel Prints	7227317	10 rolls	Box
7.5 cm x 3.6 m	Pastel Prints	7227307	10 rolls	Box
5 cm x 3.6 m	Teddy Bears	7227305	10 rolls	Box
5 cm x 3.6 m	Dinosaurs	7227302	10 rolls	Box
5 cm x 3.6 m	Soccer Balls	7227308	10 rolls	Box
7.5 cm x 3.6 m	Soccer Balls	7227309	10 rolls	Box
5 cm x 3.6 m	Camouflage	7227322	10 rolls	Box
7.5 cm x 3.6 m	Camouflage	7227306	10 rolls	Box
5 cm x 3.6 m	Paw Prints	7227320	10 rolls	Box
7.5 cm x 3.6 m	Paw Prints	7227321	10 rolls	Box
5 cm x 3.6 m	Pink Camo	7227340	10 rolls	Box
7.5 cm x 3.6 m	Pink Camo	7227341	10 rolls	Box
5 cm x 3.6 m	Winter	7227303	10 rolls	Box
7.5 cm x 3.6 m	Winter	7227304	10 rolls	Box



DELTA-CAST® SOFT

POLYESTER SEMI-RIGID CAST TAPE



Delta-Cast® Soft is a semi-rigid polyester cast tape, which allows for greater range of movement and earlier mobility to reduce some of the problems of rigid immobilization, such as muscle atrophy. It is mainly used in pediatric applications, such as club foot treatment and serial casting, as well as customized support.

FEATURES AND BENEFITS

- Smooth finish and soft edges
- 100% radiolucent
- Excellent lamination
- Provides support while allowing a limited range of motion
- Can be removed with scissors or unwrapping
- Latex Free

Dimensions	Colour	Code	Box Contents	Selling Unit
2.5 cm x 1.8 m	White	7270802	10 rolls	Box
5 cm x 3.6 m	White	7270812	10 rolls	Box
7.5 cm x 3.6 m	White	7270822	10 rolls	Box
10 cm x 3.6 m	White	7270832	10 rolls	Box
12.5 cm x 3.6 m	White	7270842	10 rolls	Box
5 cm x 3.6 m	Blue	7270805	10 rolls	Box
7.5 cm x 3.6 m	Blue	7270815	10 rolls	Box
10 cm x 3.6 m	Blue	7270825	10 rolls	Box
5 cm x 3.6 m	Red	7270808	10 rolls	Box
7.5 cm x 3.6 m	Red	7270818	10 rolls	Box
10 cm x 3.6 m	Red	7270828	10 rolls	Box
5 cm x 3.6 m	Black	7270801	10 rolls	Box
7.5 cm x 3.6 m	Black	7270811	10 rolls	Box
10 cm x 3.6 m	Black	7270821	10 rolls	Box
5 cm x 3.6 m	Paw Prints	7270809	10 rolls	Box
7.5 cm x 3.6 m	Paw Prints	7270819	10 rolls	Box

DYNACAST® PRELUDE

FIBERGLASS SPLINTING SYSTEM



Dynacast® Prelude is the most cost-effective splinting system of choice for every phase of immobilization. Its easy-stretch padding covers all edges, while its roll form eliminates any wasted product.

FEATURES AND BENEFITS

- No wrong side application
- Minimum water needed to activate
- Easy stretch padding covers all edges
- Can be opened for customizing, trimming or reinforcing
- Interlocking Single Layer Technology
- Additional sizes 2.5 cm and 15 cm
- Latex Free

SUGGESTION

- Co-Plus®, Co-Plus® LF or Easifix® Cohesive bandages are ideal to hold Dynacast® Prelude splints in place and reduce the risk of slippage.



Dimensions	Code	Box Contents	Case Contents	Selling Unit
Roll Format				
2.5 cm x 4.6 m	7141800	1 roll	2 Boxes	Case
5 cm x 4.6 m	7142000	1 roll	2 Boxes	Case
7.5 cm x 4.6 m	7142100	1 roll	2 Boxes	Case
10 cm x 4.6 m	7141500	1 roll	2 Boxes	Case
12.5 cm x 4.6 m	7141600	1 roll	2 Boxes	Case
15 cm x 4.6 m	7141700	1 roll	2 Boxes	Case
Pre-Cut Splints				
7.5 cm x 30 cm	4722200	5 splints	-	Box
7.5 cm x 40 cm	4722300	5 splints	-	Box
10 cm x 30 cm	4722400	5 splints	-	Box
10 cm x 40 cm	4722500	5 splints	-	Box
10 cm x 76 cm	4722600	4 splints	-	Box
12.5 cm x 65 cm	4722700	4 splints	-	Box
12.5 cm x 76 cm	4723500	4 splints	-	Box
Dynacast® Prelude Solo - Unpadded Pre-Cut Splints				
5 cm x 25 cm	7296300	10 splints	-	Box
7.5 cm x 30 cm	7296301	10 splints	-	Box
7.5 cm x 87 cm	7296302	10 splints	-	Box
10 cm x 38 cm	7296303	10 splints	-	Box
10 cm x 75 cm	7296304	10 splints	-	Box
12.5 cm x 75 cm	7296305	10 splints	-	Box
15 cm x 30 cm	7296306	10 splints	-	Box

DELTA-CAST® PRELUDE

POLYESTER SPLINTING SYSTEM



Delta-Cast® Prelude is a new splinting system with polyester substrate. In addition to all the great advantages as Dynacast® Prelude, Delta-Cast® Prelude offers a highly comfortable padding that allows skin to breathe, as well as smooth edges.

FEATURES AND BENEFITS

- Highly comfortable water resistant padding with 3D knit construction to allow skin to breath
- Soft edges of substrate to prevent skin irritations or abrasions
- Helps to increase patient compliance
- Latex Free

Dimensions	Code	Box Contents	Selling Unit
2.5 cm x 4.6 m	7301801	1 roll	Box
5 cm x 4.6 m	7301802	1 roll	Box
7.5 cm x 4.6 m	7301803	1 roll	Box
10 cm x 4.6 m	7301804	1 roll	Box
12.5 cm x 4.6 m	7301805	1 roll	Box



GYPSONA®

PLASTER OF PARIS BANDAGES & SLABS



The traditional rich and creamy texture of Gypsona® allows intimate molding and creates a smooth, appealing finished cast, making Gypsona® the world's leading brand.

FEATURES AND BENEFITS

- Leno gauze for extra stability
- Creamy texture
- Low plaster loss
- Central plastic core
- Latex Free

Dimensions	Code	Box Contents	Case Contents	Selling Unit
Bandages - Extra Fast Setting / 2 minutes				
5 cm x 2.7 m	LPL23	24 rolls	3 boxes	Case
7.5 cm x 2.7 m	LPL33	12 rolls	6 boxes	Case
10 cm x 2.7 m	LPL43	12 rolls	6 boxes	Case
15 cm x 2.7 m	LPL63	12 rolls	6 boxes	Case
20 cm x 2.7 m	LPL83	12 rolls	4 boxes	Case
10 cm x 4.5 m	LPL45	12 rolls	6 boxes	Case
15 cm x 4.5 m	LPL65	12 rolls	4 boxes	Case
20 cm x 4.5 m	LPL85	12 rolls	4 boxes	Case
90 cm x 5.4 m	7198737	3 rolls	-	Box
Slabs - Extra Fast Setting / 2 minutes				
7.5 cm x 75 cm	L5330	50 splints	6 boxes	Case
10 cm x 75 cm	L5430	50 splints	6 boxes	Case
12.5 cm x 75 cm	L5530	50 splints	6 boxes	Case
5 ply - 10 cm x 18 m	LP542	1 roll	-	Box
5 ply - 15 cm x 18 m	LP562	1 roll	-	Box
Bandages - Fast Setting / 4 minutes				
5 cm x 2.7 m	4000	24 rolls	3 boxes	Case
7.5 cm x 2.7 m	4010	12 rolls	6 boxes	Case
10 cm x 2.7 m	4040	12 rolls	6 boxes	Case
15 cm x 2.7 m	4045	12 rolls	6 boxes	Case
10 cm x 4.5 m	4020	12 rolls	6 boxes	Case
15 cm x 4.5 m	4030	12 rolls	4 boxes	Case
20 cm x 4.5 m	4050	12 rolls	4 boxes	Case
Slabs - Fast Setting / 4 minutes				
7.5 cm x 75 cm	4060	50 splints	6 boxes	Case
10 cm x 75 cm	4070	50 splints	6 boxes	Case
12.5 cm x 75 cm	4080	50 splints	6 boxes	Case

GYPSONA® ROLL

PLASTER OF PARIS SPLINTING SYSTEM



The foam and cotton flannel design, combined with extra-fast setting plaster makes Gypsona® Roll a clinician's product of choice for low-cost splint applications.

FEATURES AND BENEFITS

- All-in-one splinting system
- Extra-fast setting
- Easy to cut for custom fitting
- Latex Free

Dimensions	Code	Box Contents	Selling Unit
10 Layers			
5 cm x 7.3 m	5701	1 roll	Box
7.5 cm x 7.3 m	5702	1 roll	Box
10 cm x 7.3 m	5703	1 roll	Box
12.5 cm x 7.3 m	5704	1 roll	Box
15 cm x 7.3 m	5705	1 roll	Box
15 Layers			
10 cm x 7.3 m	5751	1 roll	Box
12.5 cm x 7.3 m	5752	1 roll	Box
15 cm x 7.3 m	5753	1 roll	Box

OCL®

PLASTER OF PARIS BANDAGES, SLABS & SPLINTING SYSTEM



OCL® plaster of Paris non leno bandages and slabs are designed for rapid wetting and smooth finish with a 2-4 minutes set time (no central core).

OCL® Splint Roll is fast setting plaster of Paris with foam and flannel covering, for patient comfort. (Latex Free)

Dimensions	Code	Box Contents	Case Contents	Selling Unit
5 cm x 2.7 m	41923X	12 rolls	6 boxes	Case
7.5 cm x 2.7 m	41933X	12 rolls	6 boxes	Case
7.5 cm x 4.6 m	41935X	12 rolls	6 boxes	Case
10 cm x 2.7 m	41943X	12 rolls	6 boxes	Case
10 cm x 4.6 m	41945X	12 rolls	6 boxes	Case
12.5 cm x 4.6 m	41955X	12 rolls	4 boxes	Case
15 cm x 2.7 m	41963X	12 rolls	4 boxes	Case
15 cm x 4.6 m	41965X	12 rolls	4 boxes	Case
20 cm x 4.6 m	41985X	12 rolls	3 boxes	Case
30 cm x 4.6 m	419125X	12 rolls	-	Box
45 cm x 4.6 m	419185X	12 rolls	-	Box
80 cm x 4.6 m	419325X	12 rolls	-	Box
12.5 cm x 76 cm	43530X	50 splints	6 boxes	Case
Splinting system				
5 cm x 6.1 m	SR215	1 roll	-	Box
7.5 cm x 6.1 m	SR315	1 roll	-	Box
10 cm x 6.1 m	SR415	1 roll	-	Box
12.5 cm x 6.1 m	SR515	1 roll	-	Box
15 cm x 6.1 m	SR615	1 roll	-	Box

DELTA TERRY-NET™ S

SYNTHETIC TERRY CLOTH STOCKINETTE



Delta Terry-Net™ S synthetic terry cloth stockinette is constructed of non-absorbent synthetic fibers that do not hold moisture against the skin. The terry-net design provides additional padding and stretches to five times its normal width to conform to complex anatomies.

FEATURES AND BENEFITS

- Soft synthetic fibers
- Stretches five times its width
- Available in pre-cut version with built-in thumb liner
- Latex Free

Dimensions	Code	Box Contents	Case Contents	Selling Unit
Roll Format				
5 cm x 13.5 m	7230600	1 roll	2 boxes	Case
7.5 cm x 13.5 m	7230601	1 roll	2 boxes	Case
10 cm x 13.5 m	7230602	1 roll	2 boxes	Case
12.5 cm x 13.5 m	7230603	1 roll	-	Box
15 cm x 13.5 m	7230604	1 roll	-	Box
Pre-Cut Format				
5 cm x 41 cm - Short	7230802	10 stockinettes	-	Box
7.5 cm x 41 cm - Short	7230800	10 stockinettes	-	Box
5 cm x 48 cm - Long	7230803	10 stockinettes	-	Box
7.5 cm x 48 cm - Long	7230801	10 stockinettes	-	Box

DELTA-NET® / PROTOUCH®

SYNTHETIC STOCKINETTE



Delta-Net® / Protouch® synthetic stockinettes stretch to more than twice their widths to conform comfortably to complex anatomies.

FEATURES AND BENEFITS

- Synthetic fibers
- Good conformability
- Latex Free

Dimensions	Colour	Code	Box Contents	Case Contents	Selling Unit
Delta-Net®					
5 cm x 22.75 m	White	7230700	1 roll	2 boxes	Case
7.5 cm x 22.75 m	White	7230701	1 roll	2 boxes	Case
10 cm x 22.75 m	White	7230702	1 roll	2 boxes	Case
15 cm x 22.75 m	White	7230703	1 roll	-	Box
20 cm x 22.75 m	White	7230704	1 roll	-	Box
2.5 cm x 22.75 m	Black	7272300	1 roll	2 boxes	Case
5 cm x 22.75 m	Black	7272301	1 roll	2 boxes	Case
7.5 cm x 22.75 m	Black	7272302	1 roll	2 boxes	Case
10 cm x 22.75 m	Black	7272303	1 roll	2 boxes	Case
15 cm x 22.75 m	Black	7272304	1 roll	-	Box
Protouch®					
2.5 cm x 22.75 m	White	301001	1 roll	-	Box
5 cm x 22.75 m	White	301002	1 roll	-	Box
7.5 cm x 22.75 m	White	301003	1 roll	-	Box
10 cm x 22.75 m	White	301004	1 roll	-	Box
15 cm x 22.75 m	White	301006	1 roll	-	Box
20 cm x 22.75 m	White	301008	1 roll	-	Box

PROTOUCH®

COTTON STOCKINETTE



Protouch® cotton stockinette is designed for patient protection under plaster of Paris and synthetic casting materials.

FEATURES AND BENEFITS

- High-grade cotton
- Good conformability
- Latex Free

Dimensions	Code	Box Contents	Selling Unit
2.5 cm x 22.75 m	307001	1 roll	Box
5 cm x 22.75 m	307002	1 roll	Box
7.5 cm x 22.75 m	307003	1 roll	Box
10 cm x 22.75 m	307004	1 roll	Box
15 cm x 22.75 m	307006	1 roll	Box
20 cm x 22.75 m	307008	1 roll	Box

DELTA-DRY®

WATER RESISTANT CAST PADDING AND STOCKINETTE



Delta-Dry®'s unique patented technology outperforms other water resistant orthopedic products, making it the best choice for clinicians.

FEATURES AND BENEFITS

- Allows patients to shower, swim, bathe and wash hands
- Reduces odor under the cast
- Excellent conformability to body contours
- Best long term water resistance
- Provides comfortable, fast-drying casts
- Leaves the skin healthy under the cast
- Simple and quick application
- Shows resistance to cast saw blades
- Latex Free

Dimensions	Code	Box Contents	Selling Unit
Padding			
5 cm x 2.4 m	7344300	12 rolls	Box
7.5 cm x 2.4 m	7344301	12 rolls	Box
10 cm x 2.4 m	7344302	12 rolls	Box
Delta-Dry® Stockinette			
2.5 cm x 10 m	7456400	1 roll	Box
5 cm x 10 m	7456401	1 roll	Box
7.5 cm x 10 m	7456402	1 roll	Box
10 cm x 10 m	7456403	1 roll	Box
12.5 cm x 10 m	7456404	1 roll	Box
Delta-Dry® Softliner – Padded Stockinette			
2.5 cm x 10 m	7277901	1 roll	Box
5 cm x 10 m	7277902	1 roll	Box
7.5 cm x 10 m	7277903	1 roll	Box
10 cm x 10 m	7277904	1 roll	Box
12.5 cm x 10 m	7277905	1 roll	Box



SPECIALIST® 100 / SPECIALIST® 100 STERILE

COTTON CAST PADDING



Specialist® 100 has the product features cotton users desire; soft material that feathers and self-bonds to create a smooth undercast surface. Specialist® 100 is hand-tearable and applies easily.

FEATURES AND BENEFITS

- Highly absorbent
- 100% cotton
- Excellent cohesion
- Hand-tearable
- Available in sterile version
- Latex Free

Dimensions	Code	Bag/Box Contents	Case Contents	Selling Unit
Specialist® 100				
5 cm x 3.7 m	7230200	24 rolls	6 bags	Case
7.5 cm x 3.7 m	7230201	12 rolls	6 bags	Case
10 cm x 3.7 m	7230202	12 rolls	6 bags	Case
15 cm x 3.7 m	7230203	6 rolls	6 bags	Case
Specialist® 100 Sterile				
5 cm x 3.7 m	7234600	50 rolls	-	Box
7.5 cm x 3.7 m	7234601	50 rolls	-	Box
10 cm x 3.7 m	7234602	25 rolls	-	Box
15 cm x 3.7 m	7234603	25 rolls	-	Box

SOF-ROL® / SOF-ROL® STERILE

NATURAL CAST PADDING



Sof-Rol® is made of highly absorbent surgical grade rayon suitable for use around wounds. This soft natural material creates a comfortable undercast surface.

FEATURES AND BENEFITS

- Highly absorbent
- Easy application
- 100% surgical grade rayon
- Available in sterile version
- Latex Free

Dimensions	Code	Bag/Box Contents	Case Contents	Selling Unit
Sof-Rol®				
5 cm x 3.6 m	9052	24 rolls	6 bags	Case
7.5 cm x 3.6 m	9033	12 rolls	6 bags	Case
10 cm x 3.6 m	9034	12 rolls	6 bags	Case
15 cm x 3.6 m	9036	6 rolls	6 bags	Case
Sof-Rol® Sterile				
5 cm x 3.6 m	7234510	20 rolls	-	Box
7.5 cm x 3.6 m	7234511	20 rolls	-	Box
10 cm x 3.6 m	7234512	20 rolls	-	Box
15 cm x 3.6 m	7234513	20 rolls	-	Box

PROTOUCH® PLUS

SYNTHETIC CAST PADDING



Protouch® Plus is made of non-absorbent synthetic material that does not hold moisture against the skin. Its high loft and fleece density provide extra cushioning and patient protection.

FEATURES AND BENEFITS

- Excellent cushioning
- Soft conformability
- Excellent feathering
- Odour control agent
- Autoclavable
- Latex Free

Dimensions	Code	Bag Contents	Case Contents	Selling Unit
5 cm x 2.7 m	7201402	12 rolls	24 bags	Case
7.5 cm x 2.7 m	7201403	12 rolls	16 bags	Case
10 cm x 2.7 m	7201400	12 rolls	12 bags	Case
15 cm x 2.7 m	7201401	12 rolls	8 bags	Case

DELTA-ROL®

SYNTHETIC CAST PADDING



Delta-Rol® is made of non-absorbent synthetic material that does not hold moisture against the skin. Its conformable stretch allows narrow widths fitting around small anatomies without cutting and tearing.

FEATURES AND BENEFITS

- Conformable stretch
- Comfortable
- Non-absorbent synthetic fibers
- Latex Free

Dimensions	Code	Bag Contents	Case Contents	Selling Unit
5 cm x 3.6 m	6882	24 rolls	3 bags	Case
7.5 cm x 3.6 m	6883	12 rolls	3 bags	Case
10 cm x 3.6 m	6884	12 rolls	3 bags	Case
15 cm x 3.6 m	6886	6 rolls	3 bags	Case

AMERICAN ORTHOPAEDIC™

CAST REMOVAL SYSTEM



FEATURES AND BENEFITS - CAST SAW

- Specially designed for rugged use: the perfect tool with power and endurance
- Internally lubricated for longer service
- Switch is designed and placed so that it will not break if dropped
- Hex drive provides 6 blade rotations extending blade life by 50% compared to similar competitive models
- One year limited warranty



FEATURES AND BENEFITS - CAST SAW VACUUM

- Lightweight and mobile; easy movement from room to room or patient to patient
- Unique lift-out disposable filter cartridge makes disposal quicker, easier, safer and cleaner than messy bags
- Quiet operation for patient comfort
- Function that allows the saw and vacuum to be turned on simultaneously with contact of the cast

Description	Code
Cast cutter	0295200
Cast dust vacuum	0295400
Disposable filter cartridge for 0295400	TCFC1
Portable cast dust vacuum	0295300
Disposable filter cartridge for 0295300	0295175

* Note: Both vacuums may also be adapted to the following cast cutters; Stryker #840, Stryker #940, and Zimmer #8909-01. Accessories required.




FEATURES AND BENEFITS - PORTABLE CAST DUST VACUUM

- Unique disposable filter cartridge
- Self-contained unit conveniently stores cast spreader, scissors, replacement blades and the filter storage caps
- Lightweight construction and convenient handle provide easy mobility
- Versatile; attaches easily to most cast cutters

CAST SAW BLADES

FEATURES AND BENEFITS

- “Universal Hub” fits most popular brands of cast cutters, both pin and hex drive
- Hexagon drive provides six blade rotations for hexagon drive cast cutters

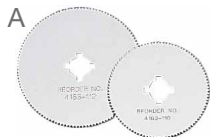

- Stainless steel
- Plaster, wood, bone



- Stainless steel
- Coated with Teflon® to reduce friction
- Plaster, synthetic cast material



- Dicronite®, unique friction reducing compound, to improve the cutting process and extend the life of the blade
- Plaster, synthetic, composites, plastics



- Titanium Nitride “Gold” Blades
- BSN’s longest lasting cast cutter blade, with coating that resists wear
- Plaster, synthetic, composites, plastics




- Stainless steel
- Special sharpening process providing sharp, long lasting cutting edge with hardened steel durability
- Plaster, composites, plastics

Code	Diameter	Bag Contents	Selling Unit
4183110 (A)	5 cm	6 blades	Bag
4183112 (A)	6.35 cm	6 blades	Bag
4183114 (B)	5 cm	1 blade	Unit
310265 (B)	6.35 cm	1 blade	Unit
0295225 (C)	5 cm	1 blade	Unit
0295226 (C)	6.35 cm	1 blade	Unit
310168 (D)	6.35 cm	1 blade	Unit
0295220 (E)	6.35 cm	1 blade	Unit

Teflon® is a registered trademark of E.I. du Pont de Nemours and Company
 Dicronite® is a registered trademark of Lubrication Sciences, Inc.


INSTRUMENTS

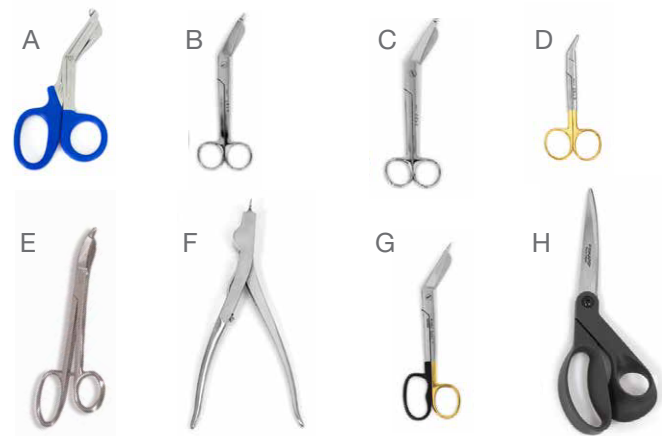

Casting instruments from BSN Medical are safe for the patient and effective for cutting casting material. They are handcrafted with high quality stainless steel to provide precision performance.

Description	Code	Selling Unit
Large utility scissors for cloth (A)	7204612	Unit
Lister scissors 180 mm (B)	7293200	Unit
Lister scissors 200 mm (C)	7293300	Unit
Small precision scissors 120 mm (D)	7293400	Unit
Bohler serrated edge scissors 235 mm (E)	7293500	Unit
Cast spreader 230 mm (F)	210136	Unit
Clean Cut™ scissors (G)	28230	Unit
Dynacast® Prelude scissors (H)	9446	Unit

This convenient cast aid provides additional patient protection during cast removal. The flexible stick slides under the cast and provides a smooth path for scissors or cast saw. Standard and metric rule adds functionality.

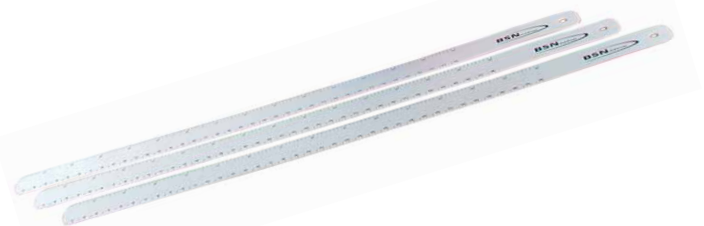

Description	Dimensions	Code	Selling Unit
Zip stick	2.5 cm x 48 cm	7204629	Unit

DELTA TERRY-NET™

ADHESIVE FLEECE EDGER



Delta Terry-Net™ adhesive fleece edger's unique stretch allows the user to pad easily around corners and contours, ensuring smooth cast edges. It also features high quality adhesive for secure placement.

FEATURES AND BENEFITS

- Highly conformable
- Unique stretch
- Secure adhesive
- Latex Free

Dimensions	Code	Colour	Bag Contents	Selling Unit
3.5 cm x 13.7 m	56012	White	1 roll	Box
3.5 cm x 13.7 m	58012	Dark Blue	1 roll	Box

DELTA TERRY-NET™

ADHESIVE PADDING



Delta Terry-Net™ adhesive padding provides added protection and cushioning for patient comfort. The adhesive backing makes attaching the padding simple and easy.

FEATURES AND BENEFITS

- Adhesive backing
- Available in terry cloth/foam and felt versions
- Latex Free

Dimensions	Code	Box Contents	Selling Unit
Terry Cloth / Foam			
58 cm x 100 cm	55012	1 sheet	Box
Felt			
58 cm x 100 cm	54012	1 sheet	Box



DELTA TERRY-NET™

ADHESIVE HOOK & STRETCH LOOP CLOSURES



This easy-to-use hook and loop system features high quality adhesive to ensure secure placement.

FEATURES AND BENEFITS

- Soft stretch loop
- High quality adhesive
- Latex Free

Dimensions	Code	Box Contents	Selling Unit
Hook			
2.5 cm x 9 m	7229500	1 roll	Box
5 cm x 9 m	7229502	1 roll	Box
Loop			
2.5 cm x 13.7 m	7229501	1 roll	Box
5 cm x 13.7 m	7229503	1 roll	Box

CASTING ACCESSORIES



This sturdy foam knee rest is used to aid in positioning a patient's limb for casting and to provide additional patient comfort during the casting process. It is covered in washable vinyl for easy clean-up.

Description	Code	Selling Unit
Knee rest	7309300	Unit



This adjustable cast stand makes applying casts easier and more comfortable for the patient. The upright can be raised or lowered from 31.8 to 53.3 cm to accommodate clinician and patient. Includes reusable plastic drip pan.

Description	Code	Selling Unit
Adjustable cast stand	58050000	Unit

CASTING ACCESSORIES



Adhesive skin traction system applied directly on the skin for non-invasive fracture reduction. Simple and effective, Tensoplast® Skin Traction Kit is used with weights (not included) and offers a reliable traction. (Latex Free)

Description	Code	Box Contents	Case Contents	Selling Unit
Skin traction kit (Adult)	7154400	1 kit	12 boxes	Case
Skin traction kit (Child)	7155400	1 kit	12 boxes	Case

These double squares walking heels with deep tread provide good patient comfort and stability by allowing regular heel-toe walking. The centre slot permits easy and secure application. (Latex Free)



Description	Dimensions	Code	Box Contents	Selling Unit
Walking heel	Adult	4183130	12 heels	Box
Walking heel	Child	4183132	12 heels	Box



Durable canvas cast shoes, with rocker bottom non-skid sole to facilitate walking, and three reliable straps. (Latex Free)

Description	Dimensions	Colour	Code	Selling Unit
Cast shoe	Extra Small	Blue	7204606	Unit
Cast shoe	Small	Blue	7204002	Unit
Cast shoe	Medium	Blue	7204001	Unit
Cast shoe	Large	Blue	7204000	Unit
Cast shoe	Extra Large	Blue	7204605	Unit
Cast shoe	Pediatric	Blue	7204003	Unit
Cast shoe with plastic loops	Extra Small	Black	7444800	Unit
Cast shoe with plastic loops	Small	Black	7444801	Unit
Cast shoe with plastic loops	Medium	Black	7444802	Unit
Cast shoe with plastic loops	Large	Black	7444803	Unit
Cast shoe with plastic loops	Extra Large	Black	7444804	Unit
Cast shoe with plastic loops	Pediatric	Black	7444600	Unit

ACTIMOVE® WALKER AIR & STANDARD

ORTHOPEDIC WALKER



The excellent features of the Actimove® Walker Air and the Actimove® Walker Standard make them unique in the market. The rigid boot, the criss-cross hook and loop straps and the full length foam ensure perfect fit and comfort. Two products, one aim – great stabilization for high patient comfort.

FEATURES AND BENEFITS

- Criss-cross ankle straps
- Liner with breathing holes (Standard version)
- Integrated air pumping system for increased stabilization (Air version)
- Latex Free

Size	Shoe Size Men	Shoe Size Women	Code	Selling Unit
Regular Low				
Extra Small	1-2.5	2-3.5	7627000	Unit
Small	3-5	4-6	7627001	Unit
Medium	5.5-9.5	6.5-10.5	7627002	Unit
Large	10-13.5	11-14.5	7627003	Unit
Extra Large	14-15	15	7627004	Unit
Regular High				
Extra Small	1-2.5	2-3.5	7627005	Unit
Small	3-5	4-6	7627006	Unit
Medium	5.5-9.5	6.5-10.5	7627007	Unit
Large	10-13.5	11-14.5	7627008	Unit
Extra Large	14-15	15	7627009	Unit
Air Low				
Extra Small	1-2.5	2-3.5	7627200	Unit
Small	3-5	4-6	7627201	Unit
Medium	5.5-9.5	6.5-10.5	7627202	Unit
Large	10-13.5	11-14.5	7627203	Unit
Extra Large	14-15	15	7627204	Unit
Air High				
Extra Small	1-2.5	2-3.5	7627205	Unit
Small	3-5	4-6	7627206	Unit
Medium	5.5-9.5	6.5-10.5	7627207	Unit
Large	10-13.5	11-14.5	7627208	Unit
Extra Large	14-15	15	7627209	Unit

ACTIMOVE® SLING

VERSATILE UPPER LIMB SUPPORT



The Actimove® Sling on a roll is a convenient support system for the upper extremity. The combination of a strap with the strong Y-tab hook closures offers users and patients an easy-to-adjust custom fit for every patient with one single product.

FEATURES AND BENEFITS

- Highly versatile; one product, many applications
- Economical; replaces several products and sizes
- Stretch resistant foam keeps arm in desired position
- Strong Y-tab hook fasteners allow easy positioning adjustments of the arm
- Strap width remains constant during use which reduces tension and increases patient comfort
- Latex Free

Description	Dimensions	Code	Box Contents	Selling Unit
Dispenser roll*	5.5 cm x 12 m	7285918	2 rolls	Box
Single pack sling **	5.5 cm x 1.9 m	7285924	1 sling	Box
Y-Tabs	-	7346600	14 hooks	Bag

* Includes 28 Y-Tab hook fasteners

** Includes 2 Y-Tab hook fasteners

CHOOSING THE RIGHT BANDAGE OR TAPE

IMMOBILIZATION

RIGID STRAPPING TAPES

- Strappal®
- Leuko® Sportstape Premium
- Leukotape® P

REHABILITATION

ELASTIC ADHESIVE TAPES AND BANDAGES

- Tensoplast® / Tensosport®
- Leukotape® K

FUNCTIONAL SUPPORT

ELASTIC BANDAGES

- Lightplast® Pro
- Co-Plus®
- Econo-San®
- Elastocrepe®
- Elastolex®

STRAPPAL®

SUPERIOR STRAPPING TAPE



Constructed from viscose fabric combined with a hypoallergenic zinc oxide adhesive, Strappal® is a high quality strapping tape offering extra strong support and efficient stabilization of joint structures.

INDICATIONS

- Rigid support for joints or strapping
- Limitation of movement
- Fixation of tubes, catheters and wound dressings

FEATURES AND BENEFITS

- Hand-tearable
- Rigid backing material
- Serrated edges
- Latex Free

Dimensions	Code	Box Contents	Selling Unit
1.25 cm x 10 m	7149615	36 rolls	Box
2.5 cm x 10 m	7148700	36 rolls	Box
4 cm x 10 m	7149000	24 rolls	Box
5 cm x 10 m	7149300	18 rolls	Box
7.5 cm x 10 m	7149401	6 rolls	Box



LEUKOTAPE® P

HIGH ADHESIVE RIGID STRAPPING TAPE



Leukotape® P is a high strength rigid strapping tape, with a strong zinc oxide adhesive. It is porous and hand-tearable, yet offers extra strong support for sprains and strains.

INDICATIONS

- Sprains & strains
- Ideal for patellofemoral taping
- Posture correction

FEATURES AND BENEFITS

- Porous
- Hand-tearable
- High strength adhesive

Dimensions	Code	Box Contents	Selling Unit
3.8 cm x 13.7 m	7616800	1 roll	Box

TENSOPLAST® / TENSOPLAST® LF

ROBUST ELASTIC ADHESIVE BANDAGE



This porous elastic adhesive bandage is made from a high-quality extensible cloth, providing firm controlled compression and strong support. The reliable adhesive keeps the bandage firmly adhered even throughout lengthy treatment.

INDICATIONS

- Firm support for injured joints
- To control effusion and swelling
- Retention of wound dressings over irregular areas and requiring compression

FEATURES AND BENEFITS

- Porous adhesive that allows skin to breathe
- Firm reliable adhesive
- Fluffy edges for comfort
- Central yellow line for accurate overlap
- Latex Free (Tensoplast® LF only)

Dimensions*	Colour	Code	Box Contents	Case Contents	Selling Unit
Tensoplast®					
5 cm x 4-5 m	Cream	7205039	12 rolls	-	Box
7.5 cm x 4-5 m	Cream	7205040	12 rolls	-	Box
10 cm x 4-5 m	Cream	7205038	12 rolls	-	Box
5 cm x 4-5 m	Cream	7205037	1 roll	12 boxes	Case
7.5 cm x 4-5 m	Cream	7205043	1 roll	12 boxes	Case
10 cm x 4-5 m	Cream	7205042	1 roll	12 boxes	Case
Tensosport®					
5 cm x 4-5 m	Bleached	7205015	36 rolls	Box	Case
7.5 cm x 4-5 m	Bleached	7205016	24 rolls	Box	Case
Tensoplast® LF - Latex Free					
5 cm x 4-5 m	Cream	7623200	12 rolls	-	Box
7.5 cm x 4-5 m	Cream	7623201	12 rolls	-	Box
10 cm x 4-5 m	Cream	7623202	12 rolls	-	Box

*Stretched

LEUKOTAPE® K

ELASTIC ADHESIVE TAPE FOR PAIN RELIEF



Leukotape® K is an elastic adhesive tape intended to help in soothing the sensation of pain and stimulating mechanoreceptors in the skin, increasing neuro-proprioceptive activity. It can also be used for encouraging lymph drainage.

INDICATIONS

- Soothing sensation of pain
- Stimulating mechanoreceptors in the skin
- Encouraging proprioception & electromyographic activity
- Encouraging lymph drainage
- Posture correction

FEATURES AND BENEFITS

- Full-day training available
- High user comfort; extremely thin and skin-friendly
- Durable for up to one week's use
- Water-resistant: can be worn in shower and during sports
- Easy to remove
- Latex Free

Dimensions	Colour	Code	Box Contents	Pack Contents	Selling Unit
2.5 cm x 5 m	Beige	7297810	1 roll	5 Boxes	Pack
5 cm x 5 m	Beige	7297811	1 roll	5 Boxes	Pack
7.5 cm x 5 m	Beige	7297812	1 roll	5 Boxes	Pack
2.5 cm x 5 m	Red	7297815	1 roll	5 Boxes	Pack
5 cm x 5 m	Red	7297816	1 roll	5 Boxes	Pack
7.5 cm x 5 m	Red	7297817	1 roll	5 Boxes	Pack
2.5 cm x 5 m	Blue	7297820	1 roll	5 Boxes	Pack
5 cm x 5 m	Blue	7297821	1 roll	5 Boxes	Pack
7.5 cm x 5 m	Blue	7297822	1 roll	5 Boxes	Pack
5 cm x 5 m	Black	7297823	1 roll	5 Boxes	Pack
5 cm x 5 m	Pink	7297825	1 roll	5 Boxes	Pack



TENSOPLAST®

ROBUST ELASTIC ADHESIVE TAPE



This porous elastic adhesive tape is made from a high-quality extensible cloth, providing firm controlled compression and strong support. The reliable adhesive keeps the bandage firmly adhered even throughout lengthy treatment.

INDICATIONS

- For secure fixation or taping with maximum conformability and comfort
- Sprains & strains requiring supportive treatment
- Retention of wound dressings over irregular areas and requiring compression
- Support bandaging of thumb and finger injuries

FEATURES AND BENEFITS

- Porous adhesive that allows skin to breathe
- Fully-spread adhesive assures firm adhesion

Dimensions*	Code	Box Contents	Case Contents	Selling Unit
2.5 cm x 4.5 m	7206902	1 roll	6 boxes	Case
5 cm x 1.5 m	7205057	1 roll	6 boxes	Case
7.5 cm x 1.5 m	72050115	1 roll	6 boxes	Case
5 cm x 4.5 m	7205218	6 rolls	-	Box
7.5 cm x 4.5 m	7205219	4 rolls	-	Box

*Stretched

TENSOPLAST®

ATHLETIC ELASTIC ADHESIVE TAPE



Tensoplast® is a light elastic adhesive tape which is suitable for medium support or compression. Its adhesive properties are long lasting and provide firm, controlled support.

INDICATIONS

- Sprains & strains
- Correction of false position & posture

FEATURES AND BENEFITS

- High adhesive properties
- Provides medium support and compression
- Water repellent
- Air permeable
- Thinner substrate

Dimensions*	Code	Box Contents	Selling Unit
5 cm x 4.5 m	04412001	24 rolls	Box
7.5 cm x 4.5 m	04413001	16 rolls	Box
10 cm x 4.5 m	04414001	16 rolls	Box

*Stretched

LIGHTPLAST® PRO

LIGHT ELASTIC ADHESIVE TAPE



Lightplast® Pro is a hand-tearable light elastic adhesive tape, suitable for light support or compression. The unique fabric backing helps to conform to body movement, and make it ideal for all-purpose taping.

INDICATIONS

- Sprains & strains
- Prophylactic use against sport injuries

FEATURES AND BENEFITS

- Hand-tearable
- Excellent adhesion
- Air permeable
- Light and close fit
- Superior elasticity

Dimensions*	Code	Box Contents	Selling Unit
2.5 cm x 4.5 m	7695200	48 rolls	Box
3.75 cm x 4.5 m	7695300	32 rolls	Box
5 cm x 4.5 m	7695400	24 rolls	Box
7.5 cm x 4.5 m	7695500	16 rolls	Box
5 cm x 6.8 m	7695600	24 rolls	Box
7.5 cm x 6.8 m	7695700	16 rolls	Box

*Stretched

CO-PLUS®/ CO-PLUS® LF

ELASTIC COHESIVE BANDAGE



Co-Plus® flexible elastic cohesive bandage provides non-slip light support and compression. Because of its unique construction, Co-Plus® adheres to itself and not to skin, aiding in patient comfort. No clips or fasteners required.

INDICATIONS

- Light support and compression
- Retention of splinting system

FEATURES AND BENEFITS

- Adheres only to itself
- Non-slip support
- Breathable
- Hand-tearable (Co-Plus® LF only)
- Latex Free (Co-Plus® LF only)

Dimensions*	Colour	Code	Box Contents	Selling Unit
Co-Plus®				
2.5 cm x 6.3 m	Beige	7210006	60 rolls	Box
5 cm x 6.3 m	Beige	7210008	36 rolls	Box
7.5 cm x 6.3 m	Beige	7210009	24 rolls	Box
10 cm x 6.3 m	Beige	7210004	18 rolls	Box
15 cm x 6.3 m	Beige	7210005	12 rolls	Box
5 cm x 6.3 m	Mixed Colours**	7210204	36 rolls	Box
7.5 cm x 6.3 m	Mixed Colours**	7210203	24 rolls	Box
10 cm x 6.3 m	Mixed Colours**	7210202	18 rolls	Box
Co-Plus® LF - Latex Free				
2.5 cm x 4.5 m	Beige	7210014	60 rolls	Box
5 cm x 4.5 m	Beige	7210016	36 rolls	Box
7.5 cm x 4.5 m	Beige	7210017	24 rolls	Box
10 cm x 4.5 m	Beige	7210018	18 rolls	Box
15 cm x 4.5 m	Beige	7210019	12 rolls	Box
5 cm x 4.5 m	Mixed Colours**	7210033	36 rolls	Box
7.5 cm x 4.5 m	Mixed Colours**	7210020	24 rolls	Box
10 cm x 4.5 m	Mixed Colours**	7210021	18 rolls	Box

*Stretched

**Mixed Colours: Neon Pink, Dark Blue, Light Blue, Purple, Neon Green and Red

ECONO-SAN®

NON-ADHESIVE ELASTIC BANDAGE



This non-adhesive bandage is made of polyester and acrylic fibres interwoven with rubber threads of high tensile strength, to provide controlled and consistent pressure.

INDICATIONS

- Sprains & strains
- Edema control (post-cast removal)
- Orthopedic support
- Post-surgical compression

FEATURES AND BENEFITS

- Consistent medium stretch
- Non-fray edges
- Excellent regain
- Non-adhesive
- Autoclavable

Dimensions*	Code	Box Contents	Selling Unit
5 cm x 4-5 m	4702	12 rolls	Box
7.5 cm x 4-5 m	4703	12 rolls	Box
10 cm x 4-5 m	4704	12 rolls	Box
15 cm x 4-5 m	4706	12 rolls	Box

*Stretched

ELASTOCREPE®

NON-ADHESIVE ELASTIC BANDAGE



Elastocrepe® is made of high quality stretch cotton/rayon fabric with non-fraying edges. The special weave helps avoid excessive stretching and provides targeted compression.

INDICATIONS

- Sprains & strains
- Edema control (post-cast removal)

FEATURES AND BENEFITS

- Lightweight and conformable
- Non-fray edges
- Washable and reusable
- Minimal extensibility, moderate stretch
- Central yellow line for accurate overlap
- Autoclavable
- Latex Free (except clips)

Dimensions*	Code	Bag Contents	Selling Unit
5 cm x 4.5 m	7208247	12 rolls	Bag
7.5 cm x 4.5 m	7208249	12 rolls	Bag
10 cm x 4.5 m	7208243	12 rolls	Bag
15 cm x 4.5 m	7208267	6 rolls	Bag

*Stretched

ELASTOLEX®

SUPERIOR NON-ADHESIVE ELASTIC BANDAGE



This non-adhesive bandage is made of high-quality fibers interwoven with rubber threads of high tensile strength, to provide controlled and consistent pressure.

INDICATIONS

- Compression of heavily swollen sprains & strains
- Edema control (post-cast removal)
- Orthopedic support
- Compression dressing
- Post-surgical compression

FEATURES AND BENEFITS

- High quality cloth for longer wear
- Non adhesive
- High compression for maximum support
- Long stretch
- Autoclavable

Dimensions*	Code	Box Contents	Selling Unit
5 cm x 4-5 m	3802H	12 rolls	Box
7.5 cm x 4-5 m	3803H	12 rolls	Box
10 cm x 4-5 m	3804H	12 rolls	Box
15 cm x 4-5 m	3806H	12 rolls	Box

*Stretched

TRICOFIX®

TUBULAR MESH BANDAGE



Tricofix® protects skin surface under compression, casts or lymphedema wrapping. Its mesh construction increases breathability and conformability.

INDICATIONS

- Protection of skin under compression
- Retention of wound dressings

FEATURES AND BENEFITS

- 100% cotton for non-irritative comfort
- Highly absorbent and well tolerated
- Easy application, may be cut anywhere, does not run or fray
- Both lengthwise and crosswise stretch for quick, easy application
- Conforms to joints securely and will not wrinkle or slip
- Durable and economical, may be washed and reused
- Autoclavable
- Latex Free

Dimensions	Code	Box Contents	Case Contents	Selling Unit
6 cm x 20 m	02195	1 roll	10 boxes	Case
8 cm x 20 m	02197	1 roll	10 boxes	Case
10 cm x 20 m	02198	1 roll	10 boxes	Case
12 cm x 20 m	02199	1 roll	10 boxes	Case

ELASTOGRIP®

TUBULAR ELASTIC SUPPORT BANDAGE



Elastogrip® is a polyester and spandex tubular elastic compression bandage which provides consistent and even compression that requires no taping.

INDICATIONS

- Medium to high support and compression
- Post cast removal

FEATURES AND BENEFITS

- Continuous uniform compression (20-25 mmHg)
- Durable elasticity
- Non-slip
- Washable and reusable
- Latex Free

Dimensions	Size	Code	Box Contents	Selling Unit
15-20 cm circumference x 10 m	5	61601	1 roll	Box
18-25 cm circumference x 10 m	6	61602	1 roll	Box
23-30 cm circumference x 10 m	7.5	61603	1 roll	Box
30-40 cm circumference x 10 m	10	61604	1 roll	Box
38-50 cm circumference x 10 m	12.5	61605	1 roll	Box

TENSOGRIP®

TUBULAR ELASTIC SUPPORT BANDAGE



Tensogrip® is an elasticized, tubular support bandage, ideal for the effective management of strains, sprains and soft tissue injuries. It provides consistent pressure and support over varying body contours.

INDICATIONS

- Light support and compression
- Retention of wound dressings

FEATURES AND BENEFITS

- Continuous uniform support (3-15 mmHg)
- Durable elasticity
- Non-slip
- Washable and reusable

Dimensions	Code	Box Contents	Selling Unit
4.5 cm x 10 m, Size A	7151300	1 roll	Box
6.25 cm x 10 m, Size B	7151400	1 roll	Box
6.75 cm x 10 m, Size C	7151900	1 roll	Box
7.5 cm x 10 m, Size D	7151500	1 roll	Box
8.75 cm x 10 m, Size E	7152000	1 roll	Box
10 cm x 10 m, Size F	7151600	1 roll	Box
12 cm x 10 m, Size G	7152100	1 roll	Box
17.5 cm x 10 m, Size J	7151700	1 roll	Box
21.5 cm x 10 m, Size K	7152300	1 roll	Box
32.5 cm x 10 m, Size L	7152200	1 roll	Box
37.5 cm x 10 m, Size M	7152400	1 roll	Box

BSN MEDICAL BRANDS

Actimove® Sling	50	Easifix® / Easifix® Cohesive	10
Actimove® Walker	49	Econo-San®	56
American Orthopaedic™	44	Elastocrepe®	56
Anti-Em/GP™	31	Elastogrip®	58
Artiflex®	30	Elastolex®	57
Cast Removal Instruments	45	Elastomull®	9
Cast Saw Blades	45	Gelocast®	26
Casting Accessories	47	Gypsona®	38
CompriFoam®	30	Hypafix®	6
Comprilan®	24	Hypafix® transparent	7
Co-Plus®	55	JOBST® Compri2	25
Coverlet®	13	JOBST® Comprifore LF	25
Covermed®	14	JOBST® FarrowWrap®	28
Coverplast® Barrier	15	JOBST® foam	30
Coverplast® Classic	14	JOBST® UlcerCARE™	27
Coverplast® Standard	15	Leukomed®	11
Cover-Roll® stretch	6	Leukomed® IV	8
Cuticell® / Cuticell® Classic	16	Leukomed® Sorbact®	13
Cuticell® Contact	17	Leukomed® T	8
Cutimed® Cavity	22	Leukomed® T plus	12
Cutimed® Gel	23	Leukoplast® Sleek	9
Cutimed® Siltec	20	Leukosan® Strip	23
Cutimed® Sorbact®	18	Leukotape® K	53
Cutimed® Sorbion®	21	Leukotape® P	52
Cutimed® Total Contact Cast Kit	22	Lightplast® Pro	55
Delta-Cast® Conformable	33	OCL®	39
Delta-Cast® Elite	34	Protouch®	40
Delta-Cast® Prelude	37	Protouch® Plus	43
Delta-Cast® Prints	35	Sof-Rol®	42
Delta-Cast® Soft	35	Specialist® 100	42
Delta-Dry®	41	Strappal®	51
Delta-Lite® Plus	32	Tensogrip®	58
Delta-Net®	40	Tensoplast®	52
Delta-Rol®	43	Tensopress®	26
Delta Terry-Net™	46	Tensosport®	52
Delta Terry-Net™ S	39	Tricofix®	57
Dynacast® Prelude	36		

BSN Medical Inc.

4455 Autoroute Laval West, Suite 255
Laval (QC) Canada H7P 4W6

www.jobstcanada.ca
www.bsnmedical.ca
Tel 1-877-978-5526
Fax 1-877-978-9703

PE160816CAT R3

